


veterinary technology for life  
**Eickemeyer**<sup>®</sup>

*promo*  
smart veterinary savings

[www.eickemeyer.com](http://www.eickemeyer.com)

Telefon +49 7461 96 580 0

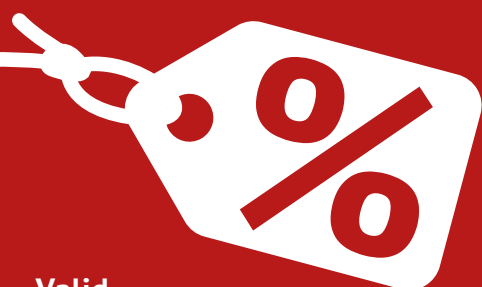

**Year End Offers  
2025**

Valid  
15/10/2025 - 31/01/2026

# EICKEMEYER® NARKOVET® II ANAESTHESIA MACHINE

## Your Veterinary Anaesthesia Solution

The EICKEMEYER® NarkoVet II is a state-of-the-art anaesthesia machine specifically designed for veterinary medicine. This compact and visually appealing device features a clear and organized structure of all components, ensuring ease of use. With the included mounting plate, an optional O<sub>2</sub> concentrator can be conveniently installed directly on the stable, mobile stand, optimising space and functionality in your practice.

- ▶ O<sub>2</sub> Flowmeter (0.1 – 4.0 l/min)
- ▶ Oxygen bypass (O<sub>2</sub> flush)
- ▶ Isoflurane vaporiser, key filled
- ▶ CO<sub>2</sub> absorber lime container with quick-change system
- ▶ Infinitely adjustable pressure relief valve
- ▶ Pressure gauge to display the system pressure
- ▶ External fresh gas outlet for open system
- ▶ Easy change of circuits using a toggle switch
- ▶ Hanger for hygienic storage of patient tubes
- ▶ Holder for anaesthetic gas filter
- ▶ Incl. O<sub>2</sub> concentrator tray, ventilator tray (NarkoVet SAV II), patient tubes (coaxial / open), rebreathing bags (0.5 l, 1 l & 2 l) and anaesthetic gas filter with scavenging hose

**213130**

## Wall mount for NarkoVet II

The wall mount allows the NarkoVet II to be mounted on the wall instead of on the mobile stand, saving space.

**21313002**

## Tray for NarkoVet II

- ▶ For monitoring devices
- ▶ Dimensions (cm): 29.5 × 19.5

**21313003**

## NarkoVet Anaesthetic Starter Kit

Complete Set, consisting of:

- ▶ Safety fill neck for isoflurane
- ▶ 3 soda lime, 4,5 kg, with color change
- ▶ Breathing bag, 3 l
- ▶ 6 plastic breathing masks (size 1 – 6)
- ▶ Complete set of EickTube tracheal tubes, 23 pcs. (2,0–10,0 mm)
- ▶ Laryngoscope ECONOMY

**213131**



## EICKEMEYER® NARKOVET SAV II SMALL ANIMAL VENTILATOR

**ORDER  
ONLINE**

[www.eickemeyer.com](http://www.eickemeyer.com)

The main cause of serious anaesthetic incidents in veterinary practices is the patient's unnoticed apnoea. Routinely controlled ventilation helps to drastically reduce the risk of anaesthesia and thus ensures the optimal safety of your patients. Ventilation is implemented via a high-performance electronic turbine, which means that no drive gas (e.g. from an external compressor) is required. The settings on the NarkoVet SAV are extremely user-friendly, just set the patient's body weight and the intelligent software automatically adjusts all other ventilation parameters.

- ▶ Electronically driven by a high-performance turbine
- ▶ Precise ventilation with a minimum tidal volume down to 4 ml
- ▶ Once inputting animal weight, smart respiratory parameters can be called up
- ▶ Ventilation modes: Volume-controlled ventilation (VCV), pressure-controlled ventilation (PCV) and Positive end-expiratory pressure (PEEP)
- ▶ Multiple alarms and clear handling prompts ensure high patient safety
- ▶ Clear "at a glance" ventilation parameters thanks to the large 7-inch colour display
- ▶ Auto system leak test with graphical guidance for quick start
- ▶ Ventilation rate: 2 – 60 bpm
- ▶ Tidal volume: 4 – 1,500 ml
- ▶ Maximum ventilation pressure: 5 – 50 cmH<sub>2</sub>O
- ▶ Release pressure -10.0 – -4.0 cmH<sub>2</sub>O
- ▶ Patient weight range: 0.4 – 150 kg
- ▶ Dimensions (in cm): L 29.1 x W 33.6 x H 46.2
- ▶ Weight: approx. 8.5 kg

**213020**



## OXYVET IV OXYGEN CONCENTRATOR

- ▶ Saving costs and effort: No oxygen bottle required, lower power consumption
- ▶ Simple operation thanks to the clear arrangement of the setting functions and monitoring information
- ▶ High patient safety through visual and acoustic alarms
- ▶ Noise level: ≤ 55 dbA
- ▶ Flow rate continuously adjustable from 0.5 – 5 l/min
- ▶ Output of up to 93 % oxygen (± 3 %)
- ▶ Dimensions (in cm): L 34 x W 37 x H 67
- ▶ Weight: 15.3 kg
- ▶ Incl. Bubble humidifier and O<sub>2</sub> hose
- ▶ 1 year warranty

**213180**



# EICKEMEYER® NARKOVET PRO ANAESTHETIC UNIT

## Maximum Safety for Your Patients

The main cause of serious anaesthesia incidents in veterinary practices is the patient's unnoticed apnoea. Routinely controlled ventilation helps to drastically minimise this risk, ensuring optimal safety for your patients. The EICKEMEYER® NarkoVet PRO features an integrated ventilator powered by an electronic turbine, eliminating the need for drive gas from an external compressor. To protect your patients and staff, the system performs an automatic leak test at start-up.

## Anaesthesia Machine

- ▶ O<sub>2</sub> Flowmeter (0.1 – 4.0 l/min)
- ▶ Oxygen bypass (O<sub>2</sub> flush)
- ▶ Isoflurane vaporiser, key filled
- ▶ CO<sub>2</sub> absorber lime container with quick-change system
- ▶ Infinitely adjustable pressure relief valve
- ▶ Pressure gauge to display the system pressure
- ▶ External fresh gas outlet for open system
- ▶ Easy change of circuits using a toggle switch
- ▶ Weighing scale for anaesthetic gas filter

## Ventilator

- ▶ Electronically controlled, no compressed air source required
- ▶ All ventilation parameters on a large, 8" touchscreen
- ▶ Patient safety through acoustic and visual alarms
- ▶ Automatic calculation of ventilation volume based on weight
- ▶ Patients from 0.5 – 150 kg can be ventilated
- ▶ Tidal volume from 5 ml, for precise ventilation of tiny patients
- ▶ Supported ventilation modes:
  - Volume-controlled ventilation (VCV)
  - Volume-supported (VS)
  - Pressure-controlled ventilation (PCV)
  - Synchronised intermittent mandatory ventilation – pressure-controlled (SIMV-PC)
  - Continuous positive airway pressure / pressure-supported ventilation (CPAP/PSV)
  - Pressure-regulated volume-controlled (PRVC)
- ▶ Synchronised intermittent mandatory ventilation – Pressure-regulated volume-controlled (SIMV-PRVC)
- ▶ Pressure-supported ventilation protection (PSVPro)
- ▶ Ventilation rate: 4 – 100 bpm
- ▶ Tidal volume: 5 – 1,500 ml
- ▶ Maximum ventilation pressure: 3 – 50 cmH<sub>2</sub>O
- ▶ Trigger pressure -20.0 to -1 cmH<sub>2</sub>O
- ▶ Mains and battery operation, battery life min. 2 hours

## Further Specs

- ▶ Incl. O<sub>2</sub> concentrator tray and patient tube with rebreathing bag (1 l)
- ▶ Dimensions (in cm): L 64 x W 67 x H 140
- ▶ Weight: approx. 43 kg
- ▶ Sevoflurane option available upon request

213040

## EICKEMEYER® NarkoVet PRO + AIR

- ▶ With second flowmeter for medical air

213041

## CO<sub>2</sub> Module for NarkoVet PRO

### The Smart Safety Upgrade for Veterinary Anaesthesia

The new CO<sub>2</sub> Module for the NarkoVet PRO system enables automated adjustment of ventilation based on the patient's exhaled CO<sub>2</sub> levels. The system continuously analyses CO<sub>2</sub> levels and automatically regulates ventilation parameters accordingly. This not only enhances monitoring precision but also ensures that respiratory support is dynamically adapted to the patient's current physiological condition – improving safety and efficiency during anaesthesia.

213045



213040

## LIFEVET CP CAPNOGRAPH AND PULSE OXIMETER

**ORDER  
ONLINE**

[www.eickemeyer.com](http://www.eickemeyer.com)

The LifeVet CP handheld Capnograph and Pulse Oximeter measures and displays the expiratory CO<sub>2</sub> concentration as well as the respiratory rate, pulse rate and oxygen saturation in dogs, cats and horses. This compact and portable device is very easy to use with its touchscreen display and intuitive interface.

- ▶ High resolution 4.3" touchscreen display
- ▶ Side-stream capnography
- ▶ EtCO<sub>2</sub> concentration: 0 – 150 mmHg
- ▶ Respiratory rate displayed in rpm
- ▶ S<sub>p</sub>O<sub>2</sub> range: 1 – 100 %
- ▶ Pulse rate: 25 – 250 bpm
- ▶ Audible and visual alarms with adjustable volume
- ▶ Long-life lithium battery: up to 18 h operating time
- ▶ Dimensions (in cm): W 7.7 x H 14.1 x D 3.6
- ▶ Weight: 420 g
- ▶ Includes: CO<sub>2</sub> tube, CO<sub>2</sub> filter, net adapter, USB cable, clip sensor

**321880**



## LIFEVET 10C MULTIPARAMETER MONITORS

The LifeVet 10C monitors are modern, powerful patient monitors that are characterised by their high level of user-friendliness, reliable technology and extensive data storage. Thanks to the large display, the vital functions can also be read easily from a distance. It is operated conveniently via touchscreen. Hot keys also provide quick access to the most important functions.

- ▶ Multi-parameter display of respiratory rate, pulse rate, oxygen saturation, blood pressure (NIBP), ECG, temperature and capnography (EtCO<sub>2</sub>) in sidestream mode
- ▶ High-resolution 10" LED touchscreen with up to 8 curves
- ▶ Freely selectable alarm limits and 360° visible alarm display
- ▶ Mains and battery operation up to 4 hours
- ▶ Dimensions (in cm): W 27 x H 22 x D 17
- ▶ Weight: ≤ 4.0 kg
- ▶ Incl. Li-Ion battery, NIBP cable with cuff, ECG cable with clamps, S<sub>p</sub>O<sub>2</sub> sensor with clips, temperature probe, CO<sub>2</sub> adapter, CO<sub>2</sub> cable and water trap

**LifeVet 10C**  
**321920**

**LifeVet 10C+**  
**321921**



321921

## INFUSOVET SP SYRINGE DRIVER / PUMP

The EICKEMEYER® InfusoVet SP is a highly accurate and reliable syringe pump. The clearly arranged, colour coded buttons make it particularly easy to use.

- ▶ 2.5" backlit LCD display
- ▶ Compatible with standard 5, 10, 20, 30, 50 and 60 ml syringes
- ▶ Safety when injecting: precise control of flow rate and injection volume
- ▶ High patient safety through acoustic alarms
- ▶ Quick start function
- ▶ With clamp for IV pole

### Technical specifications

- ▶ Mains and battery operation, battery life up to 6 hours
- ▶ Splash proof (IP34)
- ▶ Dimensions (in cm): L 33 x D 12.6 x H 11
- ▶ Weight: approx 1.8 kg

402095



## INFUSOVET IP II INFUSION PUMP

The advanced InfusoVet IP II Infusion Pump combines the latest technology with a user-friendly design to make the administration of fluids and medications more precise and efficient.

- ▶ User-friendly thanks to simple programming
- ▶ Patient safety thanks to acoustic and visual alarms
- ▶ Suitable for blood transfusions
- ▶ Stand prevents the infusion line from being disconnected
- ▶ Incl. clamp for drip stand
- ▶ Side power connection does not protrude into the cage
- ▶ Compatible with most pump specific infusion lines. May require calibration.

### Technical specifications

- ▶ Flow rate: 0.1 – 2,000 ml/h
- ▶ Volume limit: 0.1 – 9,999 ml
- ▶ Bolus rate: 0.1 – 2,000 ml/h (auto + manual)
- ▶ Occlusion level: 5 levels selectable (150 – 1,125 mmHg)
- ▶ Mains and battery operation for up to 4.5 hours
- ▶ Splash protection in accordance with IP44
- ▶ Dimensions (in cm): W 13 x H 18 x D 12
- ▶ Weight: < 1.2 kg

402087



## EICKWARM AIR WARMING SYSTEM FOR SMALL ANIMALS

**ORDER  
ONLINE**

www.eickemeyer.com

The EickWarm Air Warming System quietly delivers forced air through optionally available warming mats that are specifically designed to evenly deliver diffused warm air to the patient.

- ▶ Safe and effective management of patient body temperature
- ▶ Optimal performance through uniform heat distribution
- ▶ Extremely powerful and quiet ( $\leq 55$  dB)
- ▶ Temperature set points: Ambient air, 32 °C, 35 °C, 38 °C and 43 °C
- ▶ Three selectable fan levels (air flow)
- ▶ Reliable measurement of the hose outlet temperature
- ▶ Multiple alarm options ensure maximum patient safety
- ▶ Pole clamp for secure fixation
- ▶ Long heat hose (2 m) with flexible fastening clamp
- ▶ Compact dimensions (in cm): L 30 x W 37 x H 34
- ▶ Weight (including heat hose): 6.5 kg
- ▶ Warming mats available separately

**648000**

### Trolley for EickWarm Air

- ▶ With lockable smooth running castors and basket
- ▶ Dimensions (in cm): L 41 x W 50 x H 65

**648010**



648000

## MOECK REUSABLE WARMING MATS

These Reusable Warming Mats are placed under the patient and – thanks to their unique 3D textile structure – ensure optimal heat supply. The mats fit the body's contour perfectly, making temperature control possible in any position and without interruption at the contact points.

- ▶ Warm air is distributed evenly under the patient
- ▶ Pressure relieving positioning
- ▶ Reusable and washable at 60 °C

**Small, 55 x 65 cm**

**648110**

**Medium, 100 x 55 cm**

**648115**



648115

### MOECK Laundry Bags

We recommend using a laundry bag when washing your blankets. Warranty may be invalid if bags are not purchased and used on each wash.

**For small mat**

**648111**

**For medium / large mat**








**648121**

# EICKEMEYER® OPERATING / TREATMENT TABLES

The EICKEMEYER® X-Base Operating Table is our most popular surgery table. All components are of the highest quality, making this table a great investment in the future of your practice.

- ▶ X-Base with 4 individually adjustable feet
- ▶ Provides best possible legroom and stability
- ▶ Good working posture whether sitting or standing
- ▶ X-Aluminium slim lift column with hydraulic or mains operated pump
- ▶ Available with stainless steel X-base or heavy square base
- ▶ Self-contained, maintenance-free column
- ▶ Height adjustable from 72 – 112 cm
- ▶ Unidirectional 90° tilt mechanism, one hand control
- ▶ Stainless steel top with rounded, anti-drip edges and drain hole
- ▶ Two tabletop widths: 50 and 60 cm
- ▶ Optional heated tabletop
- ▶ Various colour options available

**602012 / 601555**

Illustration	RAL No.	Colour
		Stainless Steel
	3000	Flame Red
	6018	Lime Green
	5022	Night Blue
	2001	Red Orange
	6005	Moss Green
	4010	Pink

## Table Tie Down Cleats

- ▶ Stainless steel construction
- ▶ Mounted under the table

**602005**

## Table Mat

- ▶ For tabletop: 50 x 130 cm
- ▶ For tabletop: 60 x 130 cm

**602150**

**602160**

## Standard Rail System

Enables quick, safe and variable fixing of the patient on the table.

Consisting of:

- ▶ 2x Standard rail lengthwise
- ▶ 2x Standard rail across
- ▶ 4x Terminal block with quick fix clamps
- ▶ 10x Fixing rope, L 1 m


**605050**

## Curved Dental Tray

- ▶ With removable tray with grid
- ▶ Tray size (in cm): W 50 x D 35

**605102**

**Configure  
ONLINE**

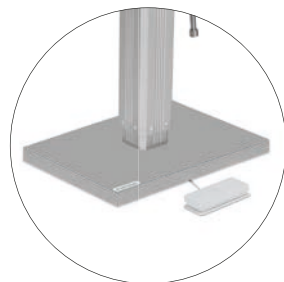

Heavy Square Base



X-Base with Dental Tray



## EICKEMEYER® SCISSOR TABLE II



The EICKEMEYER® Scissor Table II is perfect for raising and transporting larger patients without unnecessary handling and movement. The table is available in 8 different colours and the optional add-ons make this table a real all-rounder for your practice.

- ▶ Ideal for examination, prep, dental treatment and procedures
- ▶ Minimum height of 28 cm allows large dogs to walk onto the table
- ▶ Maximum height of 105 cm enables more comfortable examination of smaller patients
- ▶ All parts made from high-quality German stainless steel
- ▶ Smooth, polished tabletop (60 x 130 cm) with rounded corners, anti-drip edges and a drain hole
- ▶ Silent electric operation via foot pedal
- ▶ Includes Power Supply and Tie Down Cleats

### 604200

Illustration	RAL No.	Colour
		Stainless Steel
	3000	Flame Red
	6018	Lime Green
	5022	Night Blue
	2001	Red Orange
	6005	Moss Green
	4010	Pink
	9010	Pure White

### Additional Options

In addition to the stationary table, you can choose one or more of the following table options:

#### Battery Kit

- ▶ Includes Battery and Charger

#### Mobile Kit

- ▶ Includes Wheels, Battery and Charger

#### Integrated Scale Option

- ▶ Weighing range 0 – 120 kg
- ▶ Not compatible with Tilt Option

#### Electric Tilt Option

- ▶ Up to 35°
- ▶ Not compatible with Scale Option

#### Curved Dental Tray

- ▶ With removable tray with grid
- ▶ Tray size (in cm): W 50 x D 35

### 604303



## SCORE® AMAZONE / JUMPER STOOLS

The ever so popular SCORE® treatment stools now come in fresh colours – real eye-catchers for your practice!

- ▶ Comfortable, cushioned seat
- ▶ Easily washable synthetic leather
- ▶ Smooth running double castors
- ▶ Aluminium base diameter: 54 cm

### SCORE® JUMPER / AMAZONE Balance

The uniquely adjustable Balance mechanism from SCORE® ensures that the saddle moves subtly along with your movements. These 'micro' movements of your pelvis ensure that muscle groups and intervertebral discs are not constantly under severe static strain. This prevents pain, helps you maintain your working posture longer and exercises your core muscles. The degree of movement can be set by your own preference and adjusted by rotating the ring. The balance chairs are also available with lumbar support. This ensures optimal back support and a healthy working posture.



6101010011



6101010030

Configure  
ONLINE



Illustration	Colour
	Black
	Anthracite
	Capri Blue
	Forest Green
	Burgundy
	Orange
	Pink
	Ochre Yellow



JUMPER (wide)



AMAZONE (narrow)



Balance (flexible)



6101010045

## SCORE® MEDICAL TREATMENT AND OPERATING CHAIR

**ORDER  
ONLINE**

[www.eickemeyer.com](http://www.eickemeyer.com)

- ▶ Comfortable, padded seat with Ø 40 cm
- ▶ Easily washable synthetic leather cover
- ▶ Height adjustable by gas spring from approx. 55 – 74 cm
- ▶ Smooth PU double castors
- ▶ Optimal freedom of movement
- ▶ Aluminum base with Ø 54 cm

**610120**

Illustration	Colour
	Black
	Anthracite
	Capri Blue
	Forest Green
	Burgundy
	Orange
	Pink
	Ochre Yellow

### With lumbar support

- ▶ Compact, functional lumbar support
- ▶ Height and inclination easily adjustable
- ▶ Colour: black

**610190**

### With lumbar support and foot operation

- ▶ Foot operation allows easy adjustment of seat height during operations (operator remains sterile)
- ▶ Colour: black

**610195**



61012501

610195

**Configure  
ONLINE**



## SCORE® ERGO SUPPORT CHAIR

The SCORE® Ergo Support Chair is specially developed for ophthalmologists. It allows an ergonomic sitting position and offers you comfortable support during long periods of seated work.

Special features of the SCORE® Ergo Support chair are its ergonomically shaped armrests. They can be brought close to the body and provide your forearms optimal support. The height adjustment of the armrests allows an upright sitting position when using a microscope or other precision instrument. The Ergo Support helps you to assume a relaxed posture: the neck, shoulders and back are relieved.

- ▶ Easy height adjustment, seat and back angle
- ▶ SCORE® Ergo Shape seat with rounded corners
- ▶ Compact lumbar support can be easily positioned in the lumbar region
- ▶ Strong aluminium base with large wheels
- ▶ Height adjustment (with foot control): 52 – 71 cm
- ▶ Working height: 80 – 90 cm
- ▶ Seat dimensions: Ø 40 cm
- ▶ Foot cross dimensions: Ø 68 cm

**610060**



610060

## DR. MACH LED 130 PLUS / 130 F EXAMINATION AND THEATRE LIGHT

The DR. MACH LED 130 has been our best selling LED light over the past years with hundreds of installations across the UK and Europe. Whether in prep, dental or in the operating theatre, the LED 130 supports veterinary professionals with superior technology, invariable reliability and great design.

- ▶ Colour rendering index:  $R_a \geq 95$
- ▶ Electric dimmer: 50–100 %
- ▶ Colour temperature: 4,300 kelvin
- ▶ Light head diameter: 33 cm
- ▶ Number of LEDs: 19
- ▶ Power consumption: only 28 W
- ▶ Optional sterilisable handle (Item no. 611266)

### DR. MACH LED 130 Plus

- ▶ 100,000 lux at 1 m working distance
- ▶ Fixed light field: 12 cm at 1 m

DR. MACH LED 130 Plus	
Item No.	Model
611736	Stand Model, mobile
611737	Wall Model with 171 cm reach
611738	Ceiling Model*, up to 3 m**

### DR. MACH LED 130 F

- ▶ 70,000 lux at 1 m working distance
- ▶ Focusable light field: 14–25 cm at 1 m

DR. MACH LED 130 F	
Item No.	Model
611733	Stand Model, mobile
611734	Wall Model with 171 cm reach
611735	Ceiling Model*, up to 3 m**



611737



## DR. MACH LED 150 THEATRE LIGHT

Based on the continued success of the LED 130, the LED 150 offers a greater light head diameter and light intensity. The handle can be detached for sterilisation and additional handles can be purchased separately.

- ▶ Colour rendering index:  $R_a \geq 96$
- ▶ Electric dimmer: 50–100 %
- ▶ Light head diameter: 38 cm
- ▶ Number of LEDs: 26
- ▶ Power consumption: only 35 W

### DR. MACH LED 150 (SWING)

- ▶ 110,000 lux at 1 m working distance
- ▶ Colour temperature: 4,000 kelvin
- ▶ Fixed light field: 19 cm at 1 m

DR. MACH LED 150	
Item No.	Model
612940	Stand Model, mobile
612941	Wall Model with 171 cm reach
612942	Ceiling Model*, for ceiling heights of 2.60–3 m**

### DR. MACH LED 150 FP (SWING)

- ▶ 130,000 lux at 1 m working distance
- ▶ Colour temperature: 4,300 kelvin
- ▶ Focusable light field: 17–24 cm at 1 m

DR. MACH LED 150 FP	
Item No.	Model
612946	Stand Model, mobile
612947	Wall Model with 171 cm reach
612948	Ceiling Model*, for ceiling heights of 2.60–3 m**



612942



\* For ceiling lights, please order the ceiling ring (Item no. 612730) separately.

\*\* An additional charge applies for ceiling heights over 3 m.

## DR. MACH LED 300DF SC THEATRE LIGHT

**ORDER  
ONLINE**

[www.eickemeyer.com](http://www.eickemeyer.com)

The latest addition to the quality manufacturer DR. MACH range is the LED 300DF Theatre Light. This product offers the greatest light intensity at an extremely attractive price. The light beams of the three LED cluster bulbs can be bundled with the sterilisable handle and the light can be focused, depending on the wound size.

- ▶ 160,000 lux at 1 m working distance
- ▶ Focusable light field: 19 – 26 cm
- ▶ Colour rendering index: Ra ≥ 96
- ▶ Light head diameter: 55 cm
- ▶ Working area: 70 – 150 cm
- ▶ Number of LEDs: 36
- ▶ Lifespan of the LEDs: 60,000 h
- ▶ Electronic brightness control on the light head

### Stand Model

- ▶ On a stable  
4-foot tripod

**612926**

### Wall Model

- ▶ With a total reach  
of 175 cm

**612924**

### Ceiling Model\*

- ▶ For ceiling heights  
of 2.60 – 3 m\*\*

**612925**



612925

## DR. MACH LED 6MC (BASIC) THEATRE LIGHT

The ergonomic design of the DR. MACH LED 6MC (Basic) stands for maximum handling, minimal heat generation and optimal functionality. A highlight of the new light is the so-called cascade system. Different types of lenses ensure vertical focusing of the depth of light. This means that as the depth of the operating field increases, the light no longer needs to be readjusted manually.

- ▶ 140,000 lux at 1 m working distance
- ▶ Focusable light field: 18 – 30 cm at 1 m working distance
- ▶ Colour rendering index: Ra ≥ 98
- ▶ Number of LEDs: 69
- ▶ Lifespan of the LEDs: 60,000 h
- ▶ Electronic brightness control on the light head

### Stand Model

- ▶ On a stable  
4-foot tripod

**613010**

### Wall Model

- ▶ Including separate  
transformer box

**613011**

### Ceiling Model\*

- ▶ For ceiling heights  
up to 3 m\*\*

**613012**



613012

\* For ceiling lights, please order the ceiling ring (Item no. 612730) separately.

\*\* An additional charge applies for ceiling heights over 3 m.

## EICKTRON ELECTROSURGICAL UNIT

The EickTron is a versatile tool in veterinary medicine allowing it to be used for both monopolar and bipolar instruments. The unit's power adjustment via a foot or hand switch and its microprocessor control makes it very easy to use.

### Includes

- ▶ Handpiece with finger switch and connecting cable
- ▶ Electrode set (10 pieces)
- ▶ Neutral electrode incl. cable
- ▶ Footswitch

### Technical specifications

- ▶ Input: 220 – 240 V, 50 – 60 Hz
- ▶ Power consumption: 300 VA
- ▶ Frequency: 600 kHz
- ▶ Dimensions (in mm):  
L 288 x W 254 x H 104

**323135**



323135



323105



323106



323139

### Autoclavable Accessories

- Bipolar Forceps
- ▶ Angled
  - ▶ Length: 150 mm

**323105**

Bipolar Forceps Cable  
**323106**

Bipolar Adapter  
**323139**

## EICKTRON EVAC SMOKE EVACUATION SYSTEM

The EickTron EVAC effectively extracts and filters vapours and aerosols generated during electrosurgical procedures. It eliminates unpleasant odours, reduces bacterial and viral contamination, and ensures clear vision. The suction can be controlled automatically or manually with a footswitch. The suction power is individually adjustable.

- ▶ Maximum suction power: 1000 L/min
- ▶ Maximum suction power with filter: 300 L/min
- ▶ Maximum suction force: 250 mbar
- ▶ Filter efficiency: 99.999 %
- ▶ Dimensions (cm): W 37 x H 14.4 x 31.9
- ▶ Weight: 4 kg
- ▶ Incl. EickTron EVAC, starter filter (5-hour runtime), suction hose set, handpiece holder, power cable

**323160**

### EickTron EVAC Filter 20

The EVAC Filter 20 with ULPA and activated carbon technology offers reliable filtration of the finest particles and effective reduction of unpleasant odours and bacterial and viral contamination.

- ▶ Runtime: 20 hours
- ▶ Particle filtration: down to 0.12 micrometers
- ▶ Filter technology: ULPA with activated carbon

**32316001**



32316001



323160

## EICKVIEW KOMPAKT ENDOSKOP

**ORDER  
ONLINE**  
www.eickemeyer.com

The specially developed EickView Compact Endoscope provides the perfect solution for rhinoscopy, tracheoscopy and vaginoscopy in small animals due to its small diameter (3.8 mm) and integrated working channel.

- ▶ Outer diameter Ø 3.8 mm
- ▶ Integrated working channel Ø 1.3 mm
- ▶ Working length 14 cm
- ▶ Viewing angle: 30°
- ▶ Two integrated stopcocks
- ▶ Set includes: EickView Compact Endoscope incl. fitted otoscopy protection aid, portable LED light source incl. charger, foreign body forceps, spare battery, 2 cleaning brushes, padded carrying case

**306122**



## EICKVIEW PORTABLE LED VIDEO ENDOSCOPES

Due to their extremely compact dimensions, the EickView LED Video Endoscopes are the perfect solution for broncho-, gastro-, rhino-, tracheo-, uretero- and vaginoscopy in dogs and cats and for the examination of the upper respiratory tract in horses. The devices work independently, as the powerful LED light source and the high resolution camera of the stand-alone systems are integrated in the distal end of the endoscopes. In addition, the folding 3.5" monitor allows both direct video and image playback as well as storage via SD card. Furthermore the EickView LED Video Endoscopes can be connected to a separate monitor.

- ▶ Folding 3.5" monitor for video and picture playback
- ▶ Recording and archiving on integrated SD card
- ▶ File formats: JPG, AVI
- ▶ Integrated powerful LED light source
- ▶ Rechargeable Battery with up to 3 hours of operation
- ▶ USB port and video out
- ▶ Depth of field: 5 – 60 mm
- ▶ Direction of view: 0°
- ▶ Field of view: 110°

EickView LED Video Endoscopes					
Item No.	Description	Working length	Outer Ø	Working channel	Deflection
<b>306009</b>	EickView 60S	60 cm	3.3 mm	1.2 mm	up 180° down 130°
<b>306011</b>	EickView 60	60 cm	4.2 mm	1.3 mm	up 180° down 130°
<b>306010</b>	EickView 70	70 cm	6.0 mm	2.0 mm	up 180° down 130°
<b>306166</b>	EickView 100E	100 cm	8.0 mm	2.8 mm	up 150° down 130°
<b>306164</b>	EickView 150E	150 cm	8.0 mm	2.8 mm	up 130° down 130°



306009

# EICKVIEW 5-IN-1 ENDOSCOPY SYSTEM

## EickView 5-in-1 Endoscopy System

The EickView 5-in-1 endoscopy system is a true all-rounder for veterinary medicine. The easy-to-use, mobile device offers the option to work with both, flexible and rigid endoscopes. This means that all types of endoscopy can be covered with just one device – from diagnosis to surgical procedures.

The EickView 5-in-1 combines a large, 15.6" Full HD display, a powerful LED light source, a connection for flexible endoscopes from the EickView HD series with air and water pump as well as digital data storage via SD Card in one system.

**306200**

## Outdoor case for EickView 5-in-1

- ▶ With foam insert
- ▶ Made of robust plastic

**30620001**

## EickView HD and HD+ Video Endoscopes

The EickView HD and HD+ Video Endoscopes deliver brilliant images for bronchoscopy, gastroscopy and various other applications. The new HD+ Endoscopes are characterised by their reduced outer diameters and enlarged working channels, making them even more variable. Due to the large number of available scope sizes, they are suitable for all applications in small animals and equine.

## Technical specifications

- ▶ Viewing direction: 0°
- ▶ Distal end contains 1 CMOS Camera Chip, 2 Light Guides, Working Channel, Lens Rinsing Nozzle
- ▶ Scope of delivery: EickView Video Endoscope, Pressure Tester, 2 Cleaning Brushes, Cleaning Set, Carrying Case

## EickView 5-in-1 Camera Kit

The EickView 5-in-1 camera with zoom adapter and light guide cable completes the comprehensive package. It allows all standard rigid endoscopes to be connected to the main system.

**306201**



306201



306200

306168

EickView HD Video Endoscopes							
Art. Nr.	Description	Working length	Outer Ø	Working channel	Depth of field	Field of view	Deflection
<b>306177</b>	EickView 60 HD Endoscope	60 cm	3.5 mm	1.2 mm	3 – 50 mm	120°	Angle above: 160° Angle below: 130° Angle right / left: 0°
<b>306178</b>	EickView 100 HD Endoscope	100 cm	6.0 mm	2.0 mm	3 – 50 mm	120°	Angle above: 210° Angle below: 90° Angle right / left: 100°
<b>306266</b>	EickView 150-6 HD+ Endoscope	150 cm	6.0 mm	2.8 mm	2 – 100 mm	120°	Angle above: 210° Angle below: 90° Angle right / left: 100°
<b>306268</b>	EickView 150-8 HD+ Endoscope	150 cm	8.0 mm	2.8 mm	3 – 50 mm	120°	Angle above: 210° Angle below: 90° Angle right / left: 100°
<b>306270</b>	EickView 350 HD+ Endoscope	350 cm	13.0 mm	3.7 mm	3 – 100 mm	140°	Angle above: 180° Angle below: 180° Angle right / left: 160°



# EICKVIEW SD ENDOSCOPY SYSTEM

**ORDER  
ONLINE**

www.eickemeyer.com

## Unlock the future of veterinary care with the EickView SD Endoscopy System!

This state-of-the-art system features a high-quality endoscopy monitor with a large, intuitive touchscreen, seamlessly connecting to three versatile “semi-disposable” endoscopes. This system allows the veterinary practitioner to experience the benefits of flexible endoscopy (Gastroscopy, Bronchoscopy, Rhinoscopy, Ureteroscopy etc.) without the hefty investment of traditional video endoscopes. Each endoscope is equipped with a high-resolution CMOS image sensor and an LED light source, ensuring crystal-clear visuals.

The new EickView 65SD-URS Ureteroscope sets new standards in urology: Thanks to its small outer diameter of just 2.1 mm, the ureter can be visualised even in male dogs – a decisive advantage for minimally invasive diagnostics and treatment.

With the EickView SD endoscopy system, you can make precise diagnoses and ensure the best possible care for your veterinary patients.

### EickView SD Endoscopy Monitor

- ▶ 15.6" touchscreen with 1,920 x 1,080 Full HD resolution
- ▶ Functions: White Balance, Zoom, Freeze, photo (JPG) and video recording (MP4)
- ▶ 64 GB storage capacity
- ▶ Mains and battery operation, battery life up to 5 hours
- ▶ Outputs: HDMI, DVI, SDI and USB
- ▶ Dimensions (in cm): W 39 x H 25 x D 4
- ▶ Weight: 2.7 kg

**307000**

### EickView 150SD

- ▶ Working length: 150 cm
- ▶ Ø outside: 6.0 mm
- ▶ Ø working channel: 2.9 mm
- ▶ Direction of view: 0°
- ▶ Viewing angle: 120°
- ▶ Bending up / down: 180°
- ▶ Bending right / left: 120°

**307010**

### EickView 60SD

- ▶ Working length: 60 cm
- ▶ Ø outside: 3.2 mm
- ▶ Ø working channel: 1.2 mm
- ▶ Direction of view: 0°
- ▶ Viewing angle: 120°
- ▶ Bending up / down: 210°

**307020**

### EickView 60SD-XS

- ▶ Working length: 60 cm
- ▶ Ø outside: 2.7 mm
- ▶ Ø working channel: 1.2 mm
- ▶ Direction of view: 0°
- ▶ Viewing angle: 120°
- ▶ Bending up / down: 210°

**307030**

### EickView 65SD-URS

- ▶ Working length: 65 cm
- ▶ Ø outside: 2.1 mm
- ▶ Ø working channel: 1.2 mm
- ▶ Direction of view: 0°
- ▶ Viewing angle: 120°
- ▶ Bending up / down: 210°

**307040**

**NEW**



## MAGIC 3000 PLUS PORTABLE ULTRASOUND UNIT

### Unbeatable value for money!

The MAGIC 3000 PLUS impresses with many new functions and features. Advanced technologies including Power Doppler (PDI) help evaluate blood flow in different clinical conditions.

- ▶ High-resolution 12.1" LCD Screen, tiltable
- ▶ Cine mode with memory function (5 sec.)
- ▶ Various measurements including pregnancy calculation
- ▶ iStation™, iClear™, iZoom™, iBeam™ and iTouch™
- ▶ Second probe connection for quick probe swap
- ▶ Internal SSD storage and two USB ports
- ▶ Dimensions (in cm): W 29 x H 35 x D 17
- ▶ Lightweight: only 5.6 kg

**304045**



## PROBES FOR MAGIC 3000 PLUS / 6000 PLUS

### Micro-convex

- ▶ 5.0–8.5 MHz

**303710**

### Linear T-shaped

- ▶ 5.0–10.0 MHz

**303711**



303710



303711

### Linear Rectal

- ▶ 4.0–7.5 MHz

**304021**



304021



## MAGIC 6000 PLUS PORTABLE ULTRASOUND UNIT

### A new dimension in image quality

The MAGIC 6000 PLUS not only opens a new dimension of image quality in small animal abdominal sonography, the device can also be used in equine gynaecology and orthopaedics. With a powerful cine mode, you can always keep an overview even with restless patients under difficult conditions.

- ▶ High-resolution 15" display, tiltable
- ▶ Cine mode with memory function (5 sec.)
- ▶ PW and Power Doppler (PDI) with "Auto Trace"
- ▶ Optional Colour Doppler Upgrade (30374001)
- ▶ iStation™, iClear™, iZoom™, iBeam™, iTouch™ & iScape™
- ▶ Various application-specific measurements for abdomen and gynaecology including pregnancy calculation
- ▶ Second probe connection for quick probe swap
- ▶ Large internal SSD storage and two USB ports
- ▶ Dimensions (in cm): W 36 x H 36 x D 15
- ▶ Weight: 8.1 kg

**303740**



## CT, MRI AND ULTRASOUND POSITIONER

**ORDER  
ONLINE**

[www.eickemeyer.com](http://www.eickemeyer.com)

- ▶ Facilitates positioning of patients to allow abdominal ultrasonography
- ▶ Perfect positioning helps keep the patient calm and comfortable
- ▶ Robust, easy to clean and disinfect

### Size 1

- ▶ Suitable for cats/small dogs
- ▶ Dimensions (in cm):  
W 25 x H 10 x L 80

**303991**

### Size 3

- ▶ Suitable for large-sized dogs
- ▶ Dimensions (in cm):  
W 50 x H 20 x L 130

**303993**

### Size 2

- ▶ Suitable for medium-sized dogs
- ▶ Dimensions (in cm):  
W 35 x H 14 x L 100

**303992**



303992

## ULTRASOUND TABLE ATTACHMENT WITH CUTOUT

The redesigned EICKEMEYER® Ultrasound Table is the perfect solution for cardiac ultrasound examinations.


- ▶ Redesigned ergonomic cut out for cardiac ultrasound examinations
- ▶ 5 foldable feet for easy storage
- ▶ Ergonomically designed recess
- ▶ Provides ample elbow room for the examiner
- ▶ Compatible with Examination Tables  
(Width: 45 cm / Length: 100 cm)
- ▶ Dimensions (in cm): W 60 x H 19 x L 14

**303976**

### Ultrasound Table Mat with Cutout

- ▶ For Ultrasound Table (item no. 303976)
- ▶ Velcro straps for a secure fit

**30397602**

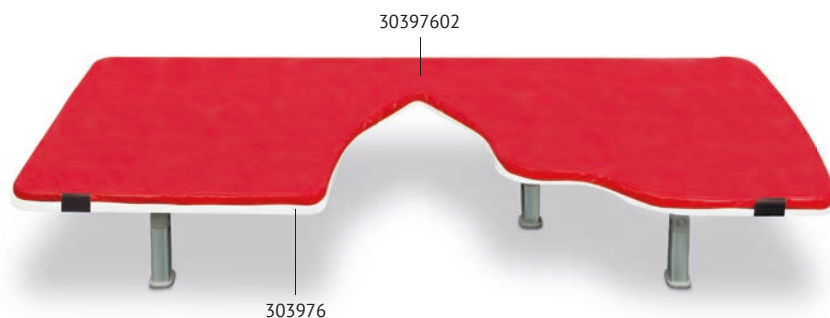

## PROFIDENT PLUS DENTAL UNIT

The Profident PLUS is the result of many years of experience. It offers a wide range of applications for both, general practitioners and dental specialists. It is characterised by its modern design combined with high mobility. The stand is adjustable in height and all handpieces are arranged within easy reach of the operator. When a handpiece is removed from its holder, the Profident PLUS automatically supplies air and water to it.

### High Speed Turbine (300,000 rpm)

- ▶ Push button burr change for friction grip burs (FG)
- ▶ Variable burr water cooling facility
- ▶ With LED Illumination

### Low Speed Air Motor (20,000 rpm) with Straight Handpiece

- ▶ For long straight shank burs (HP), diamond discs and prophy heads
- ▶ High torque, forward and reverse rotation
- ▶ Optional water cooling facility (item no. 174022)

### Built-in Dental Scaler

- ▶ Handpiece with LED Illumination
- ▶ Includes 5 scaling tips

### 3 in 1 Syringe

- ▶ Precise finger tip control of water, mist or dry air
- ▶ Water connection to connect external dental scaler (e.g. item no. 176870)

### Dental Unit Air Compressor

- ▶ Powerful, quiet, low maintenance

### Micro Oil Filter

- ▶ Helps extend the life of handpieces

**174001**

### Starter Kit

Complete set, consisting of:

- ▶ 1x Bur Kit for Contra Angle
- ▶ 1x Prophy Cups
- ▶ 1x Polishing Paste
- ▶ 1x Cleaning Brush

**EXP174004**



Scaler with LED lighting



Dental Handpiece High Speed (FG)



174001

## MOBILE DENTAL CART

### Turns every table into a dental treatment table!

With the Mobile Dental Cart from EICKEMEYER®, every operating or treatment table becomes an ergonomic dental treatment table.

- ▶ For attachment to operating and treatment tables
- ▶ Ergonomic, curved dental tray with grate
- ▶ Easily removable tray for cleaning
- ▶ Made of stainless steel
- ▶ With 4 lockable smooth running castors
- ▶ Height adjustable from 75 – 100 cm
- ▶ Tray size (in cm): W 50 x D 35

**605105**



## PROFIDENT HMS HANDPIECE MAINTENANCE SYSTEM

**ORDER  
ONLINE**  
www.eickemeyer.com

### Automated Precision for Handpiece Care

The PROFIDENT HMS is a fully automated maintenance system designed for efficient cleaning and servicing of handpieces, contra-angles, and turbines. It actively helps preserve the value of your handpieces, extends their lifespan, and ensures consistently high performance.

With its intelligent connection system, the unit accommodates up to three handpieces simultaneously. The user-friendly three-program cycle (1, 1.5, or 2 minutes) allows for customised maintenance routines tailored to your workflow.

### Technical specifications

- ▶ Maintenance capacity: up to 3 handpieces at once
- ▶ One E-type connector for low-speed handpieces and contra-angles
- ▶ Two mid west connectors for high-speed turbine handpieces
- ▶ Operates via the compressed air supply of your dental unit
- ▶ Compact dimensions (cm): W 28 × H 38 × D 27
- ▶ Weight: 7.5 kg
- ▶ Supplied with 1,000 ml maintenance oil

**174050**



## EICKEMEYER® TUB TABLE STAINLESS STEEL

This multifunctional tub table is made entirely of high-quality stainless steel and is characterised by its ergonomic design and generous storage space.

- ▶ Stylish washbasin with drawers and storage space
- ▶ Stainless steel tub with a pool depth of 15 cm
- ▶ Rounded corners for easy cleaning
- ▶ Removable, two-part grid
- ▶ Knee room for seated activities
- ▶ Faucet with extractable shower head to 70 cm
- ▶ Two water outlet options - shower and normal tap
- ▶ 3 drawers and a door
- ▶ Includes sink drainer
- ▶ Dimensions (in cm): L 140 x W 60 x H 88

**639060**



## HIRAY DENTAL X-RAY GENERATOR

---

The HiRay X-Ray Generator is the perfect solution to today's digital imaging environment and provides superb image quality.

- ▶ 0.4 mm focal spot for sharper more detailed images
- ▶ Adjustable kV (60, 65, 70)
- ▶ Adjustable time (2 sec to 10 milliseconds)

### Wall Model

- ▶ Max. reach approx. 1.98 m

**176530**

### Stand Model

- ▶ Max. reach approx. 1.3 m

**176535**



## DIGIVET CR PREMIUM DENTAL X-RAY READER

---

The new DigiVet CR PREMIUM now supports the large size 4 plate enabling rapid examination with reduced radiation exposure. The large size 4 plate can also be very useful for imaging extremities. With a resolution of 21 lp/mm, the digital imaging plate system leaves nothing to be desired.

- ▶ Optimised workflow for your practice
- ▶ Dental Imaging Plates from size 0 to 4
- ▶ Magnetic Feeding Tray with erase function
- ▶ Easy-to-use imaging Software includes DICOM (Send / Retrieve)
- ▶ Compact dimensions (in cm): H 26.8 x W 12.4 x D 35.1

Includes DigiVet CR PREMIUM Starter Kit

- ▶ 4 x Size 2 Imaging Plates (31 x 41 mm)
- ▶ 1 x Size 4c Imaging Plate (48 x 54 mm)
- ▶ 1 x Size 4 Imaging Plate (54 x 76 mm)
- ▶ 100 x Hygienic Bags, Size 2
- ▶ 100 x Hygienic Bags, Size 4c
- ▶ 100 x Hygienic Bags, Size 4

**708036**



## DENTANOMIC Set – Bronze

Complete Set, consisting of 2 handles (purple and gold) + elevator blade set (all sizes) + luxator blade set (all sizes) + sterilisation and storage rack

Up to glove size 7.5

**175050**

From glove size 8.0

**175051**

## DENTANOMIC Set – Gold

Complete Set, consisting of 5 handles (purple, gold, green, blue and pink) + elevator blade set (all sizes) + luxator blade set (all sizes) + sterilisation and storage rack + honing set + sterilisation container

Up to glove size 7.5

**175070**

From glove size 8.0

**175071**

## DENTANOMIC Set – Silver

Complete Set, consisting of 3 handles (green, blue and pink) + elevator blade set (all sizes) + luxator blade set (all sizes) + sterilisation and storage rack + honing set

Up to glove size 7.5

**175060**

From glove size 8.0

**175061**

## DENTANOMIC Set – Gold Mixed

Complete Set, consisting of 3 standard handles (green, blue and pink – up to glove size 7.5) and 2 long handles (purple and orange – from glove size 8.0) + elevator blade set (all sizes) + luxator blade set (all sizes) + sterilisation and storage rack + honing set + sterilisation container

**175072**




175070

## DENTANOMIC PERIOSTEAL ELEVATOR SET

DENTANOMIC dental instruments are an innovation in the field of dental surgery and significantly simplify tooth extraction. With the interchangeable blades, the instruments can be quickly and easily adapted to the respective surgical situation. In addition, once the blade wears out, the entire instrument no longer needs to be thrown away. The blades are extra hardened and have a special, very sharp grind.

- ▶ Set consisting of 4 periosteal elevator blades (2, 4, 6 and 8 mm) and sterilization container

**175115**

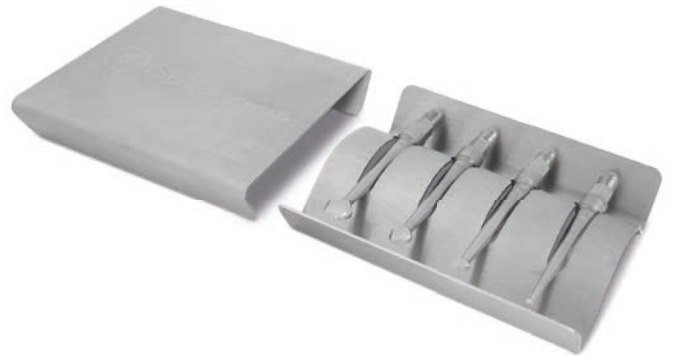

## DENTANOMIC® EQUINE SET

Following their remarkable success in the small animal market, DENTANOMIC® is proud to unveil the Full Equine Set. Developed in collaboration with Dr. Travis Henry (USA), these cutting-edge dental instruments bring DENTANOMIC®'s innovative technology to equine practitioners around the world.

The Set comes with a full range of 14 Equine Luxation Blades (3 & 5 mm in short, medium and long sizes, with straight, left and right contra-angle shapes). DENTANOMIC® Equine blade edges are narrower than the periodontal ligament space, making luxation much easier.

### DENTANOMIC® thin edged Equine blades can be used in:

- ▶ Molar Extension Bar: Articulated to allow sideways cutting motion and simultaneous rotation to follow tooth contours. An aid for molar extraction.
- ▶ New Equine Handles: Optimised for a vertical hold and increased grip. An aid for incisor extraction.

### DENTANOMIC® Equine Set Includes:

- ▶ 4x Equine Handles (purple, gold, green and blue)
- ▶ 1x Equine Handle Tray
- ▶ 1x Molar Extension Bar
- ▶ 14x Equine Luxation Blades
- ▶ 1x Blade Tray

### DENTANOMIC® Equine Set – Standard Handles

- ▶ Glove size ≤ 7.5

**175160**

### DENTANOMIC® Equine Set – Long Handles

- ▶ Glove size < 8

**175161**





## EICKEMEYER® PAWSITIONER

ORDER  
ONLINE

www.eickemeyer.com

The Pawsitioner is designed to secure patients for surgery and hands-free diagnostic imaging. It eliminates the need for tie-downs and makes it easy to clip and scrub patients in the prep room then transport them to surgery without having to re-position. The Pawsitioner is also ideal for hands-free radiographic positioning, including V-D, left or right lateral, V-D pelvis for OFA evaluation, rostral-caudal skull views and much more without re-positioning the patient for each radiograph. The cradle may also be used to hold patients for laparoscopy, CT, MRI, myelography and ultrasonography.

### Complete Set

- ▶ Consisting of 5 sizes and wall mount for space-saving storage

605920



QR code for the application video of the Pawsitioner

## EICKEMEYER® DOPPLER

Reliable and accurate non-invasive blood pressure monitor for small animals.

### EICKEMEYER® Doppler

- ▶ Compact, hand-held design
- ▶ High sensibility for a clear signal
- ▶ Adjustable sound volume and headphone connection
- ▶ Focused 8 MHz Probe with wide beam for better flow detection in the artery
- ▶ Probe adapted for positioning on the animals tail or paw
- ▶ Rechargeable battery powered, operating time: ± 6 hours
- ▶ Dimensions (in mm): L 175 x W 80 x H 25 mm
- ▶ Lightweight (300 g)

217210

### Sphygmomanometer

- ▶ Includes two cuffs: 5.0 and 7.5 cm

217309

### Doppler Blood Pressure Cuff, 1.9 cm

- ▶ Ideal for small cats

217307

### Doppler Blood Pressure Cuff, 3.3 cm

- ▶ Ideal for cats & small dogs

217308

### Doppler Blood Pressure Cuff, 10 cm

- ▶ Ideal for large dogs

217311



217210

## ICARE TONOVET PLUS TONOMETER

The iCare TONOVET Plus is easy to use. It features an advanced user interface with a clear, large display. Red/green positioning lights help to find the correct angle and aim at the middle of the cornea. The automatic measuring sequence lets you choose between taking single readings or taking all 6 measurements with one press of the button.

- ▶ Quick, painless and accurate procedure
- ▶ Easy to use and requires no calibration
- ▶ Assisting positioning lights
- ▶ Large and easy to read colour display
- ▶ Dimensions (in mm): L 215 x W 24–29 x D 35–95
- ▶ Weight: 230 g (incl. batteries)
- ▶ Includes: Transport Case, 100 Probes, 4x AA batteries

**1736007**

### Stand for TONOVET and TONOVET Plus

- ▶ Table and wall mountable

**1736008**

### Probes for TONOVET Plus

- ▶ 100 per pack

**1736006**



Application



## ICARE TONOVET PRO TONOMETER

The iCare TONOVET Pro tonometer sets new standards in intraocular pressure (IOP) measurement and is aimed at vets who have the highest demands on accuracy and ease of use. The innovative "Quick Measure" function enables faster and easier measurement of the IOP and significantly shortens the examination time, while the more responsive test tip and the optimised test tip control ensure improved handling and even more precise measurements. The 200-degree positioning freedom makes the TONOVET Pro tonometer an indispensable helper even in emergency situations. The ergonomic display is designed to allow the user to instantly read the measured values.

- ▶ Improved user-friendliness
- ▶ 200-degree positioning freedom (ideal for equine use)
- ▶ More responsive test tip
- ▶ Better test tip control improves accuracy
- ▶ Display facing the user
- ▶ Self-calibrating
- ▶ Wireless printing via Bluetooth
- ▶ Includes: Transport Case, 100 TONOVET Probes, 4x AA batteries

**173627**

### Stand for TONOVET Pro

- ▶ Table and wall mountable

**173626**

### Probes for TONOVET Pro

- ▶ 100 per pack

**1736006**



Application



## HEINE OMEGA 600 INDIRECT OPHTHALMOSCOPE

**ORDER  
ONLINE**

[www.eickemeyer.com](http://www.eickemeyer.com)

### Discover more. With less weight.

Lighter. Brighter. Better. HEINE kept the best and perfected the rest. Highly comfortable, the HEINE OMEGA 600 delivers a new class of technological innovation, along with the new and unique visionBOOST for significantly better diagnostic even in dense cataracts.

- ▶ Most precise diagnostic, up to the far retinal periphery
- ▶ Up to 20 % better view in media opacities with visionBOOST
- ▶ Lightest High-End Indirect Ophthalmoscope worldwide
- ▶ Long-lasting investment due to its robust design
- ▶ Easy disinfection thanks to smoothSURFACE

### Most precise diagnostic, up to the far retinal periphery

- ▶ Significantly better diagnostic due to unique visionBOOST
- ▶ Sharpest fundus image due to superior glass optics
- ▶ Brightest fundus image with true colour rendering thanks to our unique LEDHQ – now in even more natural colours
- ▶ Best 3D imaging due to our unique Stereoscopic Adjustment
- ▶ Technology in dilated and undilated pupils

### Includes

OMEGA 600 Indirect Ophthalmoscope (with battery and +2 D eyepiece), E4-USB power supply with cable, eyepiece 0 D

**173497**



## OPTOMED AURORA® IQ FUNDUS CAMERA

The portable, patented AURORA® IQ fundus camera is a high-end retinal imaging device for the most demanding of needs. The user-friendly technology and large screen for instant display of high-quality retinal images, even with non-dilated pupils, make this device a must-have for every ophthalmologist. The included ocular surface lens has adjustable illumination levels, color imaging and cobalt blue for fluorescence imaging. Choose between autofocus or manual operation as desired.

- ▶ Complete set including anterior module
- ▶ Non-mydratic application possible in a darkened room
- ▶ Fast image capture with autofocus or manual operation
- ▶ Real-time visibility of the retina
- ▶ Large screen
- ▶ Detailed examination with 50° field of view
- ▶ Ocular surface lens for color imaging and cobalt blue for fluorescence imaging
- ▶ Adjustable illumination levels
- ▶ Transfer to the digital patient file via WLAN and USB
- ▶ Low-noise application
- ▶ 2-year guarantee (excluding batteries)
- ▶ Scope of delivery: AURORA® IQ fundus camera, retinal module, anterior module, 2 batteries, power supply, cleaning cloth, 2 eyecups, charging station, USB cable, client software for image transfer, carrying case

### Technical specifications

- ▶ Field of view: 50° (total field of view: 50 x 40)
- ▶ Display TFT-LCD: 4", 800 x 480 pixels
- ▶ Operating time: approx. 2 hours with continuous use
- ▶ Camera dimensions (in mm): W 122 x H 202 x L 98
- ▶ Retinal module dimensions (in mm): W 68 x H 74 x L 160
- ▶ Anterior module dimensions (in mm): W 71 x H 80 x L 79
- ▶ Camera weight: 526 g

**173700**



## KOWA SL-19 PORTABLE SLIT LAMP






With its rounded shape and curved handle, the KOWA SL19 slit lamp is designed to be lightweight and ergonomic.

### The SL-19 has the following new features:

- ▶ New LED backlight that can be used simultaneously with the slit light
- ▶ New LED blue light illumination for fluorescence observation
- ▶ Infinitely adjustable slit length (1.5 – 12 mm)
- ▶ Longer operating time (up to 6 h at max. intensity)

### Technical specifications

- ▶ Magnification: 10x and 16x, working distance: about 8 – 10 cm
- ▶ Effective field of view: Ø 10 – 15 mm
- ▶ Slit width: Ø 0.1, Ø 0.2, Ø 0.8 mm
- ▶ Spot size: Ø 1, Ø 5, Ø 12 mm
- ▶ LED light intensity: 200 – 20,000 lux (continuously adjustable)
- ▶ Dimensions (in cm): W 19.7 x D 10.7 x H 23.8
- ▶ Weight (without batteries): 620 g

Illustration	Colour	Item No.
	White	173690
	Black	173691
	Blue	173692
	Green	173693
	Pink	173694

### Carrying Case for SL-19

- ▶ Dimensions (in cm): L 37 x W 26.1 x H 8.5

**173695**



173690

## KOWA SL-19 PLUS PORTABLE SLIT LAMP WITH CAMERA

The latest development from KOWA, the hand-held slit lamp SL-19 plus, impresses with its well-known ergonomic shape and light weight with an integrated, high-resolution camera. With a simple hand movement, videos can be recorded during the examination, which are automatically transmitted to an iPad via WiFi. Using special application software, these videos can be analysed and further processed on the iPad. High-resolution still images can also be extracted from the videos, which can be individually selected, saved and forwarded. With the new KOWA SL-19 plus handheld slit lamp, your examination results are always optimally analysed and documented.

### Technical specifications

- ▶ Magnification: 10x
- ▶ Working distance: approx. 100 mm
- ▶ Effective field of view: Ø 15 mm at 10x magnification
- ▶ Spot size: Ø 1, Ø 5, Ø 12 mm
- ▶ Gap length infinitely adjustable from 1.5 – 12 mm
- ▶ Light intensity: 200 – 20,000 lux (infinitely adjustable)
- ▶ Power supply: 3x AA batteries (DC 4.8 – 6.4 V)
- ▶ Dimensions (in cm): D 19.7 x W 10.7 x H 25.5
- ▶ Weight (excluding batteries): 730 g

**173696**

### Carrying Case for SL-19 plus

- ▶ Dimensions (in cm): L 37 x W 26.1 x H 8.5

**173697**



173696

## CRYOPEN XPS CRYOSURGERY DEVICE FOR OPHTHALMOLOGY

**ORDER  
ONLINE**

[www.eickemeyer.com](http://www.eickemeyer.com)

### Cryosurgery in a Nutshell

This powerful, ergonomic cryotherapy pen delivers a fine jet of nitrous oxide under high pressure, allowing you to work with millimeter precision. The CryoPen XPS is our most compact model and is suitable for intensive and versatile use, especially in ophthalmology.

- ▶ State-of-the-art microtechnology for controlled gas flow
- ▶ Ultra fine and precise
- ▶ Easy to load 8g N<sub>2</sub>O cartridges
- ▶ Aluminium housing with 2 micro-applicators
- ▶ For intensive yet gentle freezing performance
- ▶ Used for lesions, adenomas, distichiasis, epulis, eyelid adenomas, meibomian gland dysfunction, sarcoidosis, warts and other types of lesions, as well as benign tumours and carcinomas

### Technical data

- ▶ Constant treatment temperature: -89°C / -128°F
- ▶ Incl. 2 autoclavable applicators: blue, 1 – 3mm and white, 2 – 4mm
- ▶ Incl. 12x 8g cartridges
- ▶ Weight of device with full cartridge: 200 g

**173420**



## RETEVET™ PORTABLE ERG

The RETevet™ portable ERG offers an intuitive interface that allows veterinary ophthalmologists and researchers to perform electroretinograms (ERGs) to measure the retinal function of animals.

- ▶ Small, portable and battery powered
- ▶ Intuitive user interface
- ▶ Real-time results
- ▶ All flashes electronically controlled
- ▶ Easy USB connection to any PC
- ▶ Comprehensive reports in PDF formats
- ▶ Anaesthesia is optional (based on the patient)

### Technical specifications

- ▶ RGB LED technology
- ▶ 15 million colour combinations
- ▶ Dimensions (in cm): W 10.5 x H 15.5
- ▶ Weight: 240 g
- ▶ Includes: RETevet™ Device, Docking Station, 6 Jet Electrodes, 10 Needle Electrodes, Dust Cover, 3-Plug Electrode Lead, USB Cable, Power Adapter.

**173610**



## EICKLOXX SMALL OSTEOSYNTHESIS SYSTEM

The EickLoxx Small is a polyaxial locking kit for dogs and cats up to 15 kg.


The system distinguishes itself through the polyaxial placement of the custom-fit screws and combines the advantages of a locking plate system with polyaxial screw placement in  $\pm 15^\circ$  longitudinal and transverse pivoting.

Compared to conventional osteosynthesis systems, the use of locking systems (internal fixation) such as the EickLoxx Small facilitates faster healing of fractures and improved implant strength. It is for this reason that the removal of plates and screws is rarely indicated.

The biologically advantageous titanium plates are available in three sizes and can be cut to the desired length and bent in three planes. Specially designed tools are included to perform this. Titanium also reduces artefacts in post-operative imaging.

EickLoxx Small is a modular system. All implants and instruments are made in Germany using only the highest quality materials.

**185500**

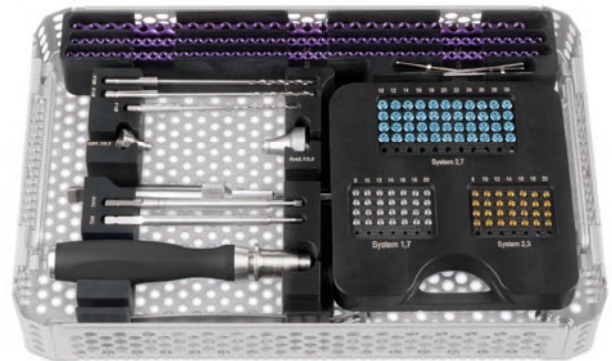

## EICKLOXX SMALL CONTOUR OSTEOSYNTHESIS SYSTEM


The EickLoxx Small Contour plate set includes various anatomically preformed osteosynthesis plates for the treatment of fractures of the distal radius and femur, ilium, acetabulum and for arthrodesis of the carpus and tarsus. It complements the cut-to-size plates of the EickLoxx Small osteosynthesis system and is therefore also designed for the treatment of cats and dogs up to 15 kg. The instruments from the EickLoxx Small osteosynthesis system are used in a similar way.

The titanium plates of the EickLoxx Small Contour plate set are available in two plate thicknesses: The 2.2 mm plates (yellow, green and magenta) can be locked polyaxially with the 1.7 mm and 2.3 mm screws and the 3.0 mm plates (light blue) optionally with the 2.7 mm, 3.5 mm and 4.0 mm screws of the EickLoxx series.

Despite contouring, the panels can be bent as required using the cabinet iron and / or the panel bending pliers with rollers and shortened using the cutting pliers.

The sieve tray with the EickLoxx Small Contour plates fits into the container Art. No. 185555 on the two sieve trays of the EickLoxx Small system.

**197900**

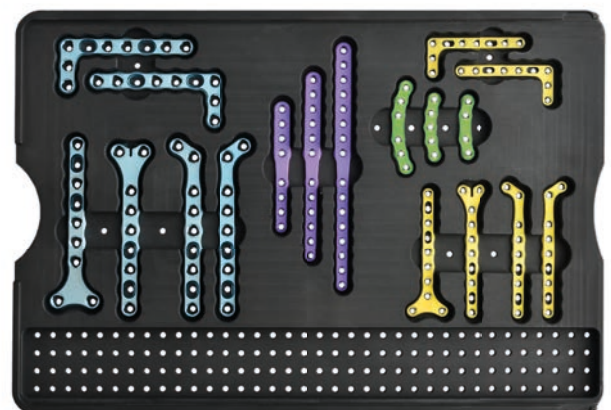

## BONE ANCHOR SET

ORDER  
ONLINE

www.eickemeyer.com

### Highest pullout strength for precise, safe ligament and tendon reconstruction

The Bone Anchor Set is an innovative, screw-thread-based anchor system for veterinary ligament and tendon reconstruction in small animals. It impresses with maximum stability, precise placement, and a sophisticated design for safe application, specifically adapted to different animal sizes, from cats to large dog breeds.

- ▶ High pullout strength: Significantly exceeds traditional hook anchors and other anchor systems
- ▶ Titanium material: Ensures high biocompatibility and infection resistance
- ▶ Self-drilling & self-tapping: Easy handling after initial perforation with a KIRSCHNER wire
- ▶ Versatile screw sizes: Diameters from 1.3 mm to 2.8 mm, available in various lengths
- ▶ Double-head drilling: Allows the insertion of multiple sutures through the same anchor head – pre- or intraoperatively using the threading aid
- ▶ Rounded head design: Minimises the risk of foreign body reactions – in contrast to conventional, classic screw heads with washers
- ▶ Ideal for collateral ligament replacement in: carpus, tarsus, elbow, knee

The complete set contains all the necessary components for immediate use – including various anchor screw sizes and a threading aid.

198200

NEW



## ORTHOVET PRO BATTERY POWERED DRILL SYSTEM

The OrthoVet PRO is a modular, battery driven drill system designed for small animal osteosynthesis. The modular system consists of a handpiece and four attachments, which can be changed intraoperatively in just a few seconds with a keyless quick coupling.

The handpiece, the Quick Release Chuck and the Jacobs Chuck are fully cannulated. The Oscillating Saw has a fivefold adjustable, variable fixation of the saw blade and allows the surgeon an optimised view of the surgical field. The AO Quick Coupling completes the system and makes the OrthoVet PRO a real all-rounder. The handpiece and attachments can be autoclaved.

- ▶ Infinitely variable speed: 0 – 1,200 rpm
- ▶ Oscillation mode: 0 – 16,000 ocs/min
- ▶ Weight (includes Battery): 933 g
- ▶ Includes: Handpiece, Jacobs Chuck, AO Chuck, Quick Release Chuck, Oscillating Saw, Saw Blades, 2 Li-ion Batteries, Charger and Transport Case

195000



## DRAGONFLY CORDLESS DRILL SYSTEM

### Compact Power for Surgical Versatility

DRAGONFLY is a modular, battery-powered drill system with four quick-release attachments for screwing, milling, drilling, and cutting. The slim handpiece features an intuitive one-handed locking mechanism. The AO attachment, quick-release chuck, and Jacobs chuck are cannulated and compatible with drill bits up to 4.0 mm. Attachments operate with variable speed, in forward/reverse and oscillation modes – up to 1,350 rpm. The system includes an oscillating saw reaching 14,000 osc/min. Handpiece and attachments are autoclavable. The central drive unit – with Li-ion battery and control electronics – is not autoclavable, but ensures long runtime and durability. DRAGONFLY delivers the performance of traditional pneumatic or wired systems – engineered and manufactured in Germany.

### Technical specifications

- ▶ Continuously adjustable speed: 0 – 1,350 rpm
- ▶ Oscillation cutting operation: 0 – 14,000 osc/min
- ▶ Battery type: Li-Ion
- ▶ Scope of delivery: DRAGONFLY handpiece including 4 attachments (AO attachment, keyless chuck for KIRSCHNER drill wires / STEINMANN nails, Jacobs drill chuck, oscillating saw), 4 saw blades with STRYKER attachment, drive unit, sterile funnel, charger and transport case

199500

**NEW**  
TPLO attachment  
available from  
the end of  
2025



## HEINE DIAGNOSTIC SET

Absolute precision is the basic requirement to ensure the creation of a high-quality product. At HEINE this principle applies to the entire production process.

Complete set, consisting of:

### NT 300 Double Charger (302550)

- ▶ Automatic charging process and overcharge protection
- ▶ Includes two NT 3.5 V Battery Handles including batteries (2x 302535)
- ▶ With short-circuit protection and a continuous ceramic brightness control

### BETA 200 LED Ophthalmoscope (300137)

The HEINE BETA 200 LED Ophthalmoscope offers a bright, clear view with durable LEDHQ technology and precise dimming, ideal for reliable mobile diagnostics.

### G100 LED Slit Illumination Head (300036)

The HEINE G100 VET LED features 1.6x magnification and precise, dimmable LED HQ lighting in a robust all-metal housing, ensuring optimal visibility and reliable diagnostics.

Replacement lamp can be ordered through item no. 300038.

### Otoscope Specula

- ▶ High quality chrome plated metal construction, which is easy to clean, disinfect and sterilise
- ▶ Metal Speculum, 4x 57 mm (300357)
- ▶ Metal Speculum, 6x 65 mm (300365)
- ▶ Metal Speculum, 7x 90 mm (300390)

EXP302550

**NT4**  
LED on  
request





## EICKVIEW OTO DIGITAL VIDEO OTOSCOPE

ORDER  
ONLINE

www.eickemeyer.com

The EickView OTO is a compact, high-performance video otoscope designed for everyday use in veterinary clinics. With its HD camera, smart WiFi connectivity, and user-friendly interface, it allows for quick and precise examinations of the ear and other small body areas – wired or wireless. The EickView OTO is designed for routine ear examinations, pre- and post-operative evaluations, as well as dermatological applications.

### Technical specifications

- ▶ HD camera: 3.9 mm lens, 6 LEDs, and 78° wide-angle view
- ▶ Magnification: 1–1.5x
- ▶ Focal Range: 10–50 mm
- ▶ Display: 1.85" IPS LCD
- ▶ Resolution: 720 x 720 px
- ▶ Built-in Storage: 8 GB
- ▶ Light adjustment: 3 levels (bright / medium / dim)
- ▶ Connection: USB-C
- ▶ Connectivity: Wireless (Android / iOS) or wired via USB-C
- ▶ Included in the Set: EickView OTO, Portable carry case, USB data & charging cable, Tongue depressor, Skin observation sleeve, Silicone ear protection sheath, Silicone sheath with 2.0 mm working channel, Biopsy forceps (1.8 mm x 280 mm)

300130



Application



## FREEZPEN 16G

### Freezpen – Cryosurgery made simple

Freezpen combines leading-edge microtechnology and ergonomic design to make precision cryosurgery easier, more precise and more accessible than ever.

- ▶ Simple, stress-free removal of skin lesions and small tumours without anaesthesia
- ▶ Simply treat lesions and tumours on the body, legs, mouth, nose, ears or in the eye area

### Technical specifications

- ▶ Treatment temperature: -89 °C / -128 °F
- ▶ Pressure: 50–55 bar
- ▶ Incl. 3 autoclavable applicators:
  - Blue: Ø 1–3 mm (230–270 seconds flow time)
  - White: Ø 2–6 mm (145–180 seconds flow time)
  - Green: Ø 4–8 mm (110–130 seconds flow time)
- ▶ Incl. 6x 16 g N<sub>2</sub>O cartridges

173400



# EOLIS AIR MANAGER

## Clean Air for Healthy Animals and People

### What is the air quality in your practice?

Air quality is often overlooked in veterinary clinics – yet it plays a crucial role in the health of patients, staff, and visitors. Whether during surgery, treatment, or routine cleaning, invisible pollutants can enter the air and pose serious health risks, potentially affecting the overall performance of your clinic.

The EOLIS Air Manager provides powerful air purification, creating a safe and hygienic environment for both animals and humans.

### Typical Airborne Contaminants in Veterinary Clinics

- ▶ VOCs from anaesthetic gases, cleaning agents, and disinfectants (e.g. isoflurane, formaldehyde)
- ▶ Pathogenic microorganisms such as bacteria, viruses, and mould spores (e.g. parvovirus, calicivirus)
- ▶ Fine and ultrafine particles from surgical procedures (e.g. electrocautery, laser)
- ▶ Odours and allergens that settle on textiles and surfaces

### Solutions with the EOLIS Air Manager

- ▶ Multi-stage filtration system: biocidal pre-filter, HEPA specialist filter, activated carbon, UV-C, and oxidation
- ▶ Deep-clean mode with active oxygen for intensive air purification in just 60 minutes
- ▶ Removes VOCs, viruses, bacteria, spores, fine dust, allergens, and odours
- ▶ Continuous air quality monitoring with real-time performance adjustment
- ▶ Ideal for treatment rooms, operating theatres, waiting areas, and staff zones
- ▶ Freestanding unit; optional wall mount or trolley available
- ▶ Made in France with a 3-year warranty

### EOLIS Air Manager 600S

- ▶ Max. coverage area (in m<sup>2</sup>): 60
- ▶ Dimensions (in cm): H 50,5 x W 38,5 x D 23
- ▶ Power consumption (min/max in W): 13 / 65
- ▶ Max. airflow (m<sup>3</sup>/h): 350
- ▶ Noise level (min/max in dB): 18 ~ 52

**503050**

### EOLIS Air Manager 1200S

- ▶ Max. coverage area (in m<sup>2</sup>): 120
- ▶ Dimensions (in cm): H 70,5 x W 39,5 x D 24
- ▶ Power consumption (min/max in W): 16 / 80
- ▶ Max. airflow (m<sup>3</sup>/h): 550
- ▶ Noise level (min/max in dB): 18 ~ 59

**503060**



503050



503060

## EICKEMEYER ICU – Precision Meets Efficiency

The new EICKEMEYER® ICU is a compact, high-end solution for advanced intensive care in veterinary practices, clinics, and ECC clinics. Two independently controlled levels enable simultaneous treatment of up to four patients – maximising capacity without increasing footprint.

Each level operates autonomously, offering precise regulation of Oxygen, CO<sub>2</sub>, Temperature, and Humidity. A digital monitoring system with smart alerts ensures maximum safety. The stainless-steel housing and integrated Ion-Based Sterilisation guarantee top-level hygiene and whisper-quiet operation – ideal for sensitive patients.

### Key Benefits

- ▶ 3D Climate Control Uniform temperature distribution creates a stable microclimate throughout the chamber.
- ▶ Real-Time O<sub>2</sub>/CO<sub>2</sub> Sensor Technology Automated gas regulation for safe, targeted therapy - compatible with central gas supply or optional oxygen concentrator.
- ▶ Negative Ion Air Purification Removes up to 98% of airborne particles, reduces pathogens and odours, promotes calm and supports immune function.
- ▶ UV Quartz Sterilisation 254 nm UV technology for effective disinfection – minimises cross-contamination.
- ▶ Therapeutic Lighting System Warm light for calming, white light for vital sign monitoring, and blue light to support cognitive function and reduce inflammation.

Whether for post-operative care, continuous oxygen therapy, or thermal management for neonates, the EICKEMEYER® ICU combines precision, hygiene, and comfort in one intelligent system. A smart investment that strengthens any modern veterinary practice.

### Technical specifications

- ▶ Touchscreen interface with real-time data
- ▶ Temperature control range: 15–38 °C
- ▶ Oxygen concentration: 21–60%
- ▶ CO<sub>2</sub> monitoring with integrated filtration
- ▶ Automatic humidity regulation
- ▶ Dimensions (cm): W 125 × D 66 × H 182
- ▶ Weight: approx. 236 kg
- ▶ Available from January 2026

**219050**

**NEW**



