veterinary technology for life

Eickemeylechnology for life

www.eickemeyer.com

Telephone +49 7461 96 580 0

Year End Offers 2024



Valid 01/10/2024 - 31/01/2025

LIFEVET PT PULSE OXIMETER WITH TEMPERATURE MEASUREMENT

The portable LifeVet PT can measure and display the functional oxygen saturation of arterial haemoglobin (S_pO_2), pulse rate and body temperature in real time in dogs, cats, and horses. All important parameter information can be easily read on the large LCD display.

- Small, handy and sturdy
- Immediately ready for use no calibration necessary
- Display of oxygen saturation from 35 to 100 %
- Display of the pulse rate from 25 to 400 beats per minute (bpm)
- Display of body temperature from 18 to 45 °C
- Alarm limits (optical & acoustic) individually adjustable
- ▶ 2.8" TFT LCD colour display with battery level indicator
- Can be connected to a computer via USB
- Test results can be stored for up to 360 hours
- ▶ Dimensions (in cm): H 14.2 x W 7.8 x D 2.8
- Weight: 250 g (device only)
- Incl. S_pO₂ sensor with two sensor clips (large & small), temperature probe, universal charger, USB cable, Li-Ion battery, protective silicone case and holder for tabletop use

321860



LIFEVET CP CAPNOGRAPH AND PULSE OXIMETER

The LifeVet CP handheld Capnograph and Pulse Oximeter measures and displays the expiratory $\rm CO_2$ concentration as well as the respiratory rate, pulse rate and oxygen saturation in dogs, cats and horses. This compact and portable device is very easy to use with its touchscreen display and intuitive interface.

- ▶ High resolution 4.3" touchscreen display
- Side-stream capnography
- EtCO₂ concentration: 0-150 mmHg
- Respiratory rate displayed in rpm
- S_pO_2 range: 1–100 %
- Pulse rate: 25 250 bpm
- Audible and visual alarms with adjustable volume
- Long-life lithium battery: up to 18 h operating time
- ▶ Dimensions (in cm): W 7.7 x H 14.1 x D 3.6
- Weight: 420 g
- ▶ Includes: CO₂ tube, CO₂ filter, net adapter, USB cable, clip sensor



LIFEVET 10M MULTIPARAMETER MONITOR

ORDER ONLINE www.eickemeyer.com

The LifeVet 10M is a fast, accurate and comfortable patient monitor that is characterised by its high ease of use and reliable technology. Its innovative technology is the ideal basis for obtaining detailed data and enabling first-class visualisations. Thanks to the large display, the vital functions can also be easily read from a distance. It is operated conveniently via touchscreen. Hot keys also provide quick access to the most important functions. The low power consumption and powerful battery make this monitor an extremely efficient solution.

- Multi-parameter display of respiratory rate, pulse rate, oxygen saturation, blood pressure (NIBP), ECG and temperature
- ▶ High-resolution 10.4" LED touchscreen with up to 7 curves
- ▶ Freely selectable alarm limits
- Mains and battery operation up to 4 hours
- Dimensions (in cm): W 32 x H 22 x D 16
- Weight: approx 3.5 kg
- Incl. Li-Ion battery, NIBP cable with cuff, ECG cable with clamps, S_pO_2 sensor with clips and temperature probe

321915



LIFEVET 10C MULTIPARAMETER MONITORS

The LifeVet 10C monitors are modern, powerful patient monitors that are characterised by their high level of user-friendliness, reliable technology and extensive data storage. Thanks to the large display, the vital functions can also be read easily from a distance. It is operated conveniently via touchscreen. Hot keys also provide quick access to the most important functions.

- Multi-parameter display of respiratory rate, pulse rate, oxygen saturation, blood pressure (NIBP), ECG, temperature and capnography (EtCO₂) in sidestream mode
- ▶ High-resolution 10" LED touchscreen with up to 8 curves
- Freely selectable alarm limits and 360° visible alarm display
- Mains and battery operation up to 4 hours
- Dimensions (in cm): W 27 x H 22 x D 17
- Weight: ≤ 4.0 kg
- Incl. Li-Ion battery, NIBP cable with cuff, ECG cable with clamps, S_pO₂ sensor with clips, temperature probe, CO₂ adapter, CO₂ cable and water trap

LifeVet 10C

LifeVet 10C+

321920 321921





EICKEMEYER® NARKOVET LIGHT ANAESTHETIC UNIT

The EICKEMEYER® NarkoVet LIGHT is a cost-effective mobile anaesthetic machine, using quality parts like the MSS Vaporiser. It includes breathing circuits and an anaesthetic gas filter.

- High accuracy Rotameter (0.1 5.0 l/min)
- EICKEMEYER® Isoflurane Vaporiser TEC 3
- Oxygen bypass (O₂ flush)
- Patient exit: Ø outside: 22 mm / Ø inside: 15 mm
- Ayre's T-Piece breathing circuit for patients under 7 kg
- Semi disposable Circle Circuit for patients over 7 kg including patient hose with Y-piece, 2.0 l rebreathing bag and APL pressure relief valve
- EickSorber anesthetic gas filter with holder and exhaust pipe
- Storage tray and hook for additional rebreathing bags
- Stable base with 5 castors (Ø 54 cm)
- Dimensions (in cm): H 142 x W 67 x D 54
- Weight: 14.5 kg

213071

OxyVet IV Oxygen Concentrator

- Saving costs and effort: No oxygen bottle required, lower power consumption
- Simple operation thanks to the clear arrangement of the setting functions and monitoring information
- High patient safety through visual and acoustic alarms
- Noise level: ≤55 dbA
- Flow rate continuously adjustable from 0.5-5 l/min
- Output of up to 93 % oxygen (± 3 %)
- Dimensions (in cm): L 34 x W 37 x H 67
- Weight: 15.3 kg
- Incl. Bubble humidifier and O₂ hose
- 1 year warranty

213180

Starter Kit for NarkoVet LIGHT

This set provides the necessary accessories to perform anaesthesia on dogs and cats.

Complete Set, consisting of:

- Isoflurane Safety Key Filler
- Soda Lime Anaesthesia, 5 l, with colour indicator
- 2 Rebreathing Bags (1 and 3 l)
- 6 Anaesthesia Masks (size 1 6)
- ▶ EickTube Tracheal Tube Set with 23 tubes (2.0 10.0 mm)
- Laryngoscope Basic







213071

EICKEMEYER® NARKOVET II ANAESTHETIC UNIT



The EICKEMEYER® NarkoVet II is a state-of-the-art anaesthesia machine specifically designed for veterinary medicine. This compact and visually appealing device features a clear and organized structure of all components, ensuring ease of use. With the included mounting plate, an optional O2 concentrator can be conveniently installed directly on the stable, mobile stand, optimising space and functionality in your practice.

- O₂ Flowmeter (0.1 4.0 l/min)
- Oxygen bypass (O₂ flush)
- Isoflurane vaporiser, key filled
- CO₂ absorber lime container with quick-change system
- Infinitely adjustable pressure relief valve
- Pressure gauge to display the system pressure
- External fresh gas outlet for open system
- Easy change of circuits using a toggle switch
- Hanger for hygienic storage of patient tubes
- Holder for anaesthetic gas filter
- Incl. O₂ concentrator tray, ventilator tray (NarkoVet SAV II), patient tubes (coaxial / open), rebreathing bags (0.5 l, 1 l & 2 l) and anaesthetic gas filter with scavenging hose

213130

OxyVet III Oxygen Concentrator (5 l)

An oxygen concentrator draws in room air, separates the oxygen from other gases in the environment and delivers the oxygen at high concentrations.

- Flowrate: 0.5 to 5 l/min
- Oxygen concentration: 93 % ±3 %
- Sound level: 40 dbA (50 Hz)
- Electrical requirements: 230 V / 50 Hz
- Power consumption: approx 290 watt @ 2 l/min and approx. 312 watt at 5 l/min
- Dimensions (in cm): L 34 x W 30.4 x H 62
- Weight: 16.3 kg
- 3 Year Warranty

213175



NarkoVet Anaesthetic Starter Kit

Complete Set, consisting of:

- Safety fill neck for isoflurane
- 3 soda lime, 4,5 kg, with color change
- Breathing bag, 3 l
- 6 plastic breathing masks (size 1 6)
- Complete set of EickTube tracheal tubes, 23 pcs. (2,0-10,0 mm)
- Laryngoscope ECONOMY









EICKEMEYER® NARKOVET PRO ANAESTHETIC UNIT

Maximum Safety for Your Patients with the EICKEMEYER® NarkoVet PRO

The main cause of serious anaesthesia incidents in veterinary practices is the patient's unnoticed apnoea. Routinely controlled ventilation helps to drastically minimise this risk, ensuring optimal safety for your patients. The EICKEMEYER® NarkoVet PRO features an integrated ventilator powered by an electronic turbine, eliminating the need for drive gas from an external compressor. To protect your patients and staff, the system performs an automatic leak test at start-up.

Anaesthesia Machine

- O₂ Flowmeter (0.1 4.0 l/min)
- Oxygen bypass (O₂ flush)
- Isoflurane vaporiser, key filled
- ▶ CO₂ absorber lime container with quick-change system
- Infinitely adjustable pressure relief valve
- Pressure gauge to display the system pressure
- External fresh gas outlet for open system
- Easy change of circuits using a toggle switch
- Holder for anaesthetic gas filter

Ventilator

- Electronically controlled, no compressed air source required
- All ventilation parameters on a large, 8" touchscreen
- Patient safety through acoustic and visual alarms
- Automatic calculation of ventilation volume based on weight
- ▶ Patients from 0.5 150 kg can be ventilated
- Tidal volume from 5 ml, for precise ventilation of the smallest patients
- Supported ventilation modes:

Volume-controlled ventilation (VCV)

Volume-supported (VS)

Pressure-controlled ventilation (PCV)

Synchronised intermittent mandatory ventilation – pressure-controlled (SIMV-PC)

Continuous positive airway pressure / pressure-supported ventilation (CPAP/PSV)

Pressure-regulated volume-controlled (PRVC)Synchronised intermittent mandatory ventilation – Pressure-regulated volume-controlled (SIMV-PRVC)

Pressure-supported ventilation protection (PSVPro)Ventilation rate: 4 - 100 bpm

- Tidal volume: 5 1,500 ml
- Maximum ventilation pressure: 3 50 cmH₂O
- Trigger pressure -20.0 to -1 cmH₂O
- Mains and battery operation, battery life min. 2 hours

Further Specs

- Incl. O_2 concentrator tray, patient tubes (coaxial / open), rebreathing bag (0.5, 1 & 2 l) and anaesthetic gas filter with scavenging hose
- ▶ Dimensions (in cm): L 64 x W 67 x H 140
- Weight: approx. 43 kg

213040

EICKEMEYER® NarkoVet PRO O2

Like NarkoVet® PRO but with a second flow meter for air 213041

NarkoVet Anaesthetic Starter Kit

Complete Set, consisting of:

- ▶ Safety fill neck for isoflurane
- 3 soda lime, 4,5 kg, with color change
- Breathing bag, 3 l
- ▶ 6 plastic breathing masks (size 1 6)
- Complete set of EickTube tracheal tubes, 23 pcs. (2,0-10,0 mm)
- Laryngoscope ECONOMY





213041

FREEZPEN

The ideal solution for cryogenic treatment of skin lesions and small tumors Powerful, precise, practical and can be used without anesthesia. This powerful, ergonomic cryotherapy pen releases a fine jet of nitrous oxide under high pressure and ensures millimeter-precision in seconds The Freezpen attachments are autoclavable or can be wiped down with a non-corrosive disinfectant solution Your advantage: Freezpen cartridges remain sealed even after use.

Freezpen 16G

- Ideal for regular use and treatment of larger or deeper lesions
- Includes $6x 16 \text{ g N}_2\text{O}$ cartridges, 1 Freezpen attachment each, blue (1-3 mm), white (2-6 mm), green (4-8 mm)

173400

Freezpen 35G

- Higher performance, for intensive use and treatment of larger or deeper lesions
- Includes $4x 35 \text{ g N}_2\text{O}$ cartridges, 1 Freezpen attachment each, green (4–8 mm), white (60 mm long, 1–3 mm) and orange (60 mm long, 2–6 mm)

173410

 N₂O Cartridges for
 N₂O Cartridges for

 Freezpen 16G
 Freezpen 35G

 ▶ Pack of 12
 ▶ Pack of 6

 173405
 173515

Attachment for Freezpen 16G

E17340001 E17340002

Attachment for Freezpen 16G and 35G

▶ Green (4-8 mm)

E17340003

Attachment for Freezpen 35G

White (60 mm long, 1-3 mm)

E17341001 E17341002

Orange (60 mm long, 2-6 mm)







173400





173410





QR code to the video of the Freezpen

INFUSOVET SP SYRINGE DRIVER / PUMP

The EICKEMEYER® InfusoVet SP is a highly accurate and reliable syringe pump. The clearly arranged, colour coded buttons make it particularly easy to use.

- 2.5" backlit LCD display
- Compatible with standard 5, 10, 20, 30, 50 and 60 ml syringes
- Safety when injecting: precise control of flow rate and injection volume
- High patient safety through acoustic alarms
- Quick start function
- With clamp for IV pole

Technical specifications

- Mains and battery operation, battery life up to 6 hours
- Splash proof (IP34)
- Dimensions (in cm): L 33 x D 12.6 x H 11
- Weight: approx 1.8 kg

402095



EICKEMEYER® I.V. FLUID WARMER

- For warming up all intravenously administered liquids
- LED display with temperature display
- Temperature control from 35 to 42 °C in 0.1 °C steps
- Warm-up time: less than 2 min
- Dimensions (in cm): L 18 x W 7 x D 3.5
- Weight: 200 g (without mounting clamp)

402092



I.V. GIVING SET

- Pressure / gravity infusion
- Individually packaged sterile sets
- DEHP & Latex free
- Compatible with infusion pump InfusoVet, InfusoVet Easy II, InfusoVet Connect and InfusoVet IP
- Please note: If you use a pump that is not from EICKEMEYER®, your pump for this infusion set must be calibrated accordingly.
- Pack of 100 pieces
- Extension 050166 is compatible with this giving set.



INFUSOVET CONNECT INFUSION PUMP



- ▶ 4.3" color touchscreen for easy operation
- Adjustable audible and visual alarms
- Suitable for bloods
- Flow rate: 0.1 2,000 ml/h
- ▶ Bolus rate: 0.1 2,000 ml/h (automatic and manual)
- Occlusion level: 5 levels selectable (50 900 mmHg)
- Mains and battery powered (up to 10 h)
- With clamp for IV pole and cage holder
- ► For EICKEMEYER® infusion sets (item no. 402071)
- Dimensions (in cm): W 13.8 x H 13.2 x D 9
- Weight: 1.5 kg

402090





INFUSOVET IP INFUSION PUMP

The EICKEMEYER® InfusoVet IP is a programmable infusion pump that combines safety, performance and reliability. The pump is easy to use and the infusion sets are easy to insert.

- 2.5" backlit LCD display
- User-friendly through simple programming
- ▶ High level of patient safety thanks to acoustic and visual alarms (in the event of air bubbles, constipation, empty battery, end of infusion)
- Also suitable for blood transfusions
- With clamp for IV pole
- ► For EICKEMEYER® infusion lines (item no. 402071)

Technical specifications

- ▶ Flow rate: 0.1 1,500 ml/h
- ▶ Volume limit: 1 9,999 ml
- ▶ Bolus rate: 0.2 600 ml/h (manual)
- Occlusion level: 3 levels selectable (150 900 mmHg)
- Mains and battery operation, battery life up to 4 hours
- Splash water protection according to IP34
- Dimensions (in cm): W 15 x H 25 x D 9
- ▶ Weight: < 1.45 kg



X-BASE OPERATING TABLE

The EICKEMEYER® X-Base Operating Table is our most popular surgery table. All components are of the highest quality, making this table a great investment in the future of your practice.

- X-Base with 4 individually adjustable feet
- Provides best possible legroom and stability
- Good working posture whether sitting or standing
- Slim lift column made of aluminium, pump hydraulic or mains operated
- Self-contained, maintenance-free lift column
- ▶ Height adjustable from 72 112 cm
- ▶ Unidirectional 90° tilt mechanism, one hand control
- Stainless steel top with rounded, anti-drip edges and drain hole
- Tabletop available in two different sizes (50 x 130 cm or 60 x 130 cm)
- Optional heated tabletop
- Various colour options available

602012

Illustration	RAL No.	Colour
		Stainless Steel
	3000	Flame Red
	6018	Lime Green
	5022	Night Blue
	2001	Red Orange
	6005	Moss Green
	4010	Pink

Table Tie Down Cleats

- Stainless steel construction
- Mounted under the table

602005

Standard Rail System

Enables quick, safe and variable fixing of the patients on the table.

Consisting of:

- 2x Standard Rails, side
- 2x Standard Rails, end
- ▶ 4x Tie Down Cleats, quick release
- ▶ 10x Tie Down Ropes, L 1 m

605050

Curved Dental Tray

- With removable tray with grid
- ▶ Tray size (in cm): W 50 x D 35





OPERATING TABLE V-TOP

The EICKEMEYER® V-Top Table features two stainless steel table panels which can be positioned independently to create the angle that holds the animal just right for all kinds of procedures. It also features an easy to clean removable gutter with a drain.

Curved Dental Tray

- With removable tray with grid
- Tub size (in cm): W 50 x D 35

602065





Technical specifications

- Adjustable top panels: 0-90°
- Stainless steel construction
- Unidirectional 75° tilt mechanism
- Adjustable height: 72 to 112 cm
- 4 Individually adjustable feet
- 3 Year Warranty
- Tabletop dimensions (in cm): L 140 x W 50

Electric-lift Mechanism 602060

Hydraulic-lift Mechanism

602070



MOBILE OPERATING TABLE X-BASE

The best selling EICKEMEYER® X-base operating and treatment tables are now available in a mobile version. Equipped with either maintenance-free, smooth-running pump hydraulics or a quiet, battery-operated electric motor with manual transmission, the table sets new standards in mobility and stability.

- Slim, attractive hydraulic lifting column made of anodized aluminium
- ▶ Height adjustable from 72 112 cm allows comfortable treatment or surgery, sitting or standing
- X-shaped substructure made of high-quality stainless steel and four smooth-running, braked double castors
- ▶ Tabletop tiltable up to 90°, easily adjustable with one hand
- With rounded corners, anti-drip edges and a drain hole
- Easy to clean and disinfect
- With smooth foot pump

50 x 130 cm 60 x 130 cm 601420 601430



LOW LIFT SCISSOR TABLE II

The EICKEMEYER® Low Lift Scissor Table II is perfect for lifting and transporting larger patients without unnecessary handling and movement. The basic stationary table is available in eight different colours and the additional options available make this table a real all-rounder for your practice.

- Ideal for examination, prep, dental treatment and procedures
- Minimum height of 28 cm allows large dogs to walk onto the table
- Maximum height of 105 cm allows you to examine smaller patients without having to hunch and lean forward
- All parts made from high quality German stainless steel
- Smooth, polished tabletop (60 x 130 cm) with rounded corners, anti-drip edges and a drain hole
- ▶ Silent electric operation via foot pedal
- Includes power supply and tie down cleats

604200

Illustration	RAL No.	Colour
		Stainless Steel
	3000	Flame Red
	6018	Lime Green
	5022	Night Blue
	2001	Red Orange
	6005	Moss Green
	4010	Pink
	9010	Pure White

Additional Options

In addition to the stationary table, you can choose one or more of the following table options:

Battery Kit

Includes Battery and Charger

Mobil Kit

Includes Wheels, Battery and Charger

Electric Tilt Option

For tabletops

Integrated Scale Option

▶ Weighing range 0–120 kg

Curved Dental Tray

- With removable tray with
- Tray size (in cm): W 50 x D 35





HYDRAULIC TABLE

ORDER ONLINE www.eickemeyer.com

The hydraulic mechanism helps protect your back. The heavy duty frame ensures stability throughout its range. All four castors can be locked simultaneously with the foot switch.

Technical specifications

- Minimum height: 49.5 cmMaximum height: 98 cmMaximum load: 100 kg
- Four lockable castors
- Tabletop dimensions (in cm): L 122 x W 60
- Includes: 6 tie-down cleats

600200



SHOR-LINE 3-KENNEL COMBINATION

- Our most popular unit with its mobile base and removable partition
- The bottom cage can be easily divided by a partition allowing 4 animals to reside comfortably
- Suitable for 3 to 4 animals
- Two Top row kennels: 610 x 610 mm
- One Bottom row kennel with double door: 1.219 x 762 mm
- With mobile substructure and removable partition
- ▶ Top row: 2 cages (each W 61 x H 61 x D 72 cm)
- Lower row: 1 cage with double door (W 121.9 x H 76.2 x D 72 cm)



SCORE® AMAZONE/JUMPER STOOLS

The ever so popular SCORE® treatment stools now come in fresh colours – real eye-catchers for your practice!

- Comfortable, cushioned seat
- Easily washable synthetic leather
- Smooth running double castors
- Aluminium base diameter: 54 cm



AMAZONE (narrow)

SCORE® JUMPER / AMAZONE Balance

The uniquely adjustable Balance mechanism from SCORE® ensures that the saddle moves subtly along with your movements. These 'micro' movements of your pelvis ensure that muscle groups and intervertebral discs are not constantly under severe static strain. This prevents pain, helps you maintain your working posture longer and exercises your core muscles. The degree of movement can be set by your own preference and adjusted by rotating the ring.



JUMPER (wide)



Balance (flexible)



Illustration	Colour
	Black
	Anthracite
	Capri Blue
	Forest Green
	Burgundy
	Orange
	Pink
	Ochre Yellow



SCORE® MEDICAL TREATMENT AND OPERATING CHAIR

- Comfortable, padded seat with Ø 40 cm
- Easily washable synthetic leather cover
- ▶ Height adjustable by gas spring from approx. 55 74 cm
- Smooth PU double castors
- Optimal freedom of movement
- Aluminum base with Ø 54 cm

610120

With lumbar support

- Compact, functional lumbar support
- ▶ Height and inclination easily adjustable
- Colour: black

610190

With lumbar support and foot operation

- Foot operation allows easy adjusting of seat height during operations (operator remains sterile)
- Colour: black

610195











610195

ZWEEZI SADDLE STOOL

- Padded synthetic leather seat
- Adjustable height: 58 85 cm
- Seat width: 40 cm
- Gas lift mechanism
- Base diameter: 54 cm
- 5 star base with swivable castors
- Colour: evergreen

610176



ZWOOZI ROUND STOOL

- Padded synthetic leather seat
- Adjustable height: 64-85 cm
- Seat width: 35 cm
- ▶ Gas lift mechanism
- Base diameter: 54 cm
- ▶ 5 star base with swivel castors
- Colour: evergreen



DR. MACH LED 130 PLUS / 130 F EXAMINATION AND THEATRE LIGHTS

The DR. MACH LED 130 has been our best selling LED light over the past years with hundreds of installations across the UK and Europe. Whether in prep, dental or in the operating theatre, the LED 130 supports veterinary professionals with superior technology, invariable reliability and great design.

Colour rendering index: R_a ≥ 95
 Electric dimmer: 50 – 100 %
 Colour temperature: 4,300 kelvin
 Light head diameter: 33 cm
 Number of LEDs: 19

Power consumption: only 28 WOptional sterilisable Handle (611266)



DR. MACH LED 130 Plus (SWING)

- ▶ 100,000 lux at 1 m working distance
- Fixed light field: 12 cm at 1 m

DR. MACH LED 130 PLUS	
Item No. Model	
611736	Stand Model, mobile
611737	Wall Model with 171 cm reach
611738	Ceiling Model*, up to 3 m

DR. MACH LED 130 F (SWING)

- 70,000 lux at 1 m working distance
- Focusable light field: 14-25 cm at 1 m

DR. MACH LED 130 F	
Item No. Model	
611733	Stand Model, mobile
611734	Wall Model with 171 cm reach
611735	Ceiling Model*, up to 3 m

^{*} For ceiling lights, please order the ceiling ring (item no. 612730) separately

DR. MACH LED 150/150 FP THEATRE LIGHTS

Based on the continuing success of the LED 130, the LED 150 offers a larger light head diameter and more light intensity. The Sterilisable Handle can be removed for sterilisation and additional handles can be purchased separately.

Colour rendering index: R_a ≥ 96
 Electric dimmer: 50 – 100 %
 Colour temperature: 4,300 kelvin
 Light head diameter: 38 cm

Number of LEDs: 26

Power consumption: only 35 W

Fixed light field: 19 cm at 1 m


DR. MACH LED 150 (SWING)

DR. MACH LED 150 FP (SWING)

- ▶ 130,000 lux at 1 m working distance
- Focusable light field: 17 24 cm at 1 m

DR. MACH LED 150	
Item No.	Model
612940	Stand Model, mobile
612941	Wall Model with 171 cm reach
612942	Ceiling Model*,
	for ceiling heights of 2.60 – 3 m

110,000 lux at 1 m working distance

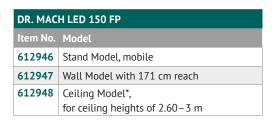



^{*} For ceiling lights, please order the ceiling ring (item no. 612730) separately

DR. MACH LED 300DF THEATRE LIGHTS

The latest addition to the range of the quality manufacturer DR. MACH impresses with the highest light intensity at an extremely attractive price. The light beams of the three LED clusters can be bundled using the sterilisable handle and, depending on the size of the wound area, also focused.

- 160,000 lux at 1 m working distance
- ▶ Focusable light field: 17 27 cm
- Colour rendering index: Ra ≥ 95
- ▶ Height adjustment: 118 cm
- Working area: 70 140 cm
- Number of LEDs: 36
- Lifespan of the LEDs: 60,000 hours
- ▶ Electronic brightness control on the luminaire body
- Including power supply

Stand Model

On a stable 4-foot tripod

612926

Wall Model

With a total reach of 175 cm

612924

Ceiling Model*

For ceiling heights of 2.60 – 3 m

612925



612925

DR. MACH LED 6MC (BASIC) THEATRE LIGHTS

The ergonomic design of the Dr. Mach LED 6sc stands for maximum handling, minimal heat generation and optimal functionality. A highlight of the new light is the so-called cascade system. Different types of lenses ensure vertical focusing of the depth of light. This means that as the depth of the operating field increases, the light no longer needs to be readjusted manually.

- ▶ 140,000 lux at 1 m working distance
- ▶ Focusable light field: 18 30 cm at 1 m working distance
- Colour rendering index: Ra ≥ 98
- Number of LEDs: 69
- Lifespan of the LEDs: 60,000 h
- ▶ Electronic brightness control on the light head

Stand Model

On a stable 4-foot tripod

613010

Wall Model

Including separate transformer box

613011

Ceiling Model*

For ceiling heights of 2.60 – 3 m

613012



^{*} For ceiling lights, please order the ceiling ring (item no. 612730) separately

^{*} For ceiling lights, please order the ceiling ring (item no. 612730) separately

MAGIC 500 PLUS PORTABLE ULTRASOUND UNIT

The MAGIC 500 PLUS sets a new standard of mobility for pregnancy testing in the field. With only 1.1 kg and an integrated rechargeable battery, it ensures maximum comfort and user friendliness.

Features

- Long service life: due to the long battery life of the battery supplied, unforeseen extended scanning sessions are not a problem
- ▶ Efficient workflow: Performing measurements on stored images ensures an efficient examination process
- Via a WiFi connection images and videos can be transmitted wirelessly to external devices (PC, laptop, tablet, smartphone) as well as monitors
- ▶ Higher diagnostic reliability with Cine Memory (256 frames)
- Individually adjustable ultrasound image
- Protection by splash-proof venting and membrane keypad
- Mounting rings and straps for easy hanging around the neck

Areas Of application

- Pregnancy examinations in cattle and horses
- Adjustment aids for optimal work
- ▶ Standard measurements and vet pregnancy calculation
- Anti-reflective ultrasound image individually adjustable
- Measurements below 1 mm possible (fat layer thickness measurement)
- Back fat and loin thickness calculation

Technical specifications

- ▶ 15 cm LED Monitor
- ▶ Black-and-white system
- Images and video clips can be saved on the device
- Incl. a switchable, changeable transducer 6.5 MHz, linear rectal, frequency switching of: 4.5 / 5.5 / 6.5 / 7.5 MHz
- Incl. a glove for the left hand (right glove can be reordered)
- Incl. Carrying case
- ▶ Integrated battery with up to 5 hours runtime
- Dimensions (in cm): W 16 x H 16 x D 7
- Weight: 1.1 kg

303790

Glove (right hand)

Right hand glove for MAGIC 500 PLUS 30379001



30379001



303790

MAGIC 3000 PLUS PORTABLE ULTRASOUND UNIT

ORDER ONLINE www.eickemeyer.com

Unbeatable value for money!

The MAGIC 3000 PLUS impresses with many new functions and features. Advanced technologies including Power Doppler (PDI) help evaluate blood flow in different clinical conditions.

- ▶ High-resolution 12.1" LCD Screen, tiltable
- Cine mode with memory function (5 sec.)
- Various measurements including pregnancy calculation
- iStation™, iClear™, iZoom™, iBeam™ and iTouch™
- Second probe connection for quick probe swap
- Internal SSD storage and two USB ports
- ▶ Dimensions (in cm): W 29 x H 35 x D 17
- Lightweight: only 5.6 kg

304045



PROBES FOR MAGIC 3000 PLUS / 6000 PLUS

Micro-convex

▶ 5.0 – 8.5 MHz

303710

Linear Rectal

▶ 4.0 – 7.5 MHz

304021

Linear T-shaped ▶ 5.0 – 10.0 MHz

303711







MAGIC 6000 PLUS PORTABLE ULTRASOUND UNIT

A new dimension in image quality

The MAGIC 6000 PLUS not only opens a new dimension of image quality in small animal abdominal sonography, the device can also be used in equine gynaecology and orthopaedics. With a powerful cine mode, you can always keep an overview even with restless patients under difficult conditions.

- High-resolution 15" display, tiltable
- Cine mode with memory function (5 sec.)
- PW and Power Doppler (PDI) with "Auto Trace"
- Optional Colour Doppler Upgrade (30374001)
- iStation™, iClear™, iZoom™, iBeam™, iTouch™ & iScape™
- Various application-specific measurements for abdomen and gynaecology including pregnancy calculation
- Second probe connection for quick probe swap
- Large internal SSD storage and two USB ports
- Dimensions (in cm): W 36 x H 36 x D 15
- Weight: 8.1 kg



CT, MRI AND ULTRASOUND POSITIONER

- Facilitates positioning of patients to allow abdominal ultrasonography
- Perfect positioning helps keep the patient calm and comfortable
- Robust, easy to clean and disinfect

Size 1

- Suitable for cats/small dogs
- Dimensions (in cm): W 25 x H 10 x L 80

303991

Size 3

- Suitable for large-sized dogs
- Dimensions (in cm): W 50 x H 20 x L 130

303993



Size 2

- Suitable for medium-sized dogs
- Dimensions (in cm): W 35 x H 14 x L 100

303992

ULTRASOUND TABLE ATTACHMENT WITH CUTOUT


The redesigned EICKEMEYER® Ultrasound Table is the perfect solution for cardiac ultrasound examinations.

- Redesigned ergonomic cut out for cardiac ultrasound examinations
- 5 foldable feet for easy storage
- Ergonomically designed recess
- Provides ample elbow room for the examiner
- Compatible with Examination Tables (Width: 45 cm / Length: 100 cm)
- Dimensions (in cm): W 60 x H 19 x L 14

303976

Ultrasound Table Mat with Cutout

- For Ultrasound Table (item no. 303976)
- Velcro straps for a secure fit

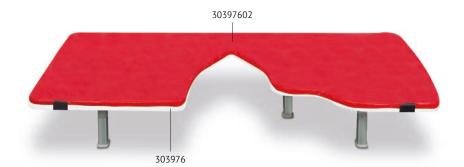

EICKEMEYER® PAWSITIONER

ORDER ONLINE www.eickemeyer.com

The Pawsitioner is designed to secure patients for surgery and hands-free diagnostic imaging. It eliminates the need for tie-downs and makes it easy to clip and scrub patients in the prep room then transport them to surgery without having to re-position. The Pawsitioner is also ideal for hands-free radiographic positioning, including V-D, left or right lateral, V-D pelvis for OFA evaluation, rostral-caudal skull views and much more without re-positioning the patient for each radiograph. The cradle may also be used to hold patients for laparoscopy, CT, MRI, myelography and ultrasonography.

Complete Set

Consisting of 5 sizes and wall mount for space-saving storage

605920



ANIMAL STRETCHER TROLLEY

We get a lot of positive feedback regarding this Animal Stretcher Trolley. It will quickly become essential for transporting patients within your practice, and you will wonder how you ever managed without it.

This mobile trolley has four lockable castors so patients can be transported and secured easily.

- 4 lockable castors
- Aluminium frame with a washable nylon cover
- Maximum load capacity: 60 kg
- ▶ Dimensions (in cm): L 120 x W 52 x H 80



HIRAY PLUS X-RAY MACHINE

One of the smallest and most powerful portable High Frequency X-Ray units in the world! The HiRay PLUS is a portable unit, which is suitable for outdoor use, as well as inside small animal practices. Special features include powerful high-frequency inverter, lightweight unit with small dimensions, and an integrated laser pointer for precise distance measuring.

Technical specifications

- ▶ High RF power: 40 to max. 100 kV / 30 mA, 2.0 kW
- kV setting: 40 to 100 kV (in 1 jV steps)
- 32 mAs settings, 0.3 to 20 mAs with illuminated LCD display
- ▶ Anatomical programming of 8 kV/mAs combinations
- Focus: 1.2 x 1.2 mm
- Collimator with bright halogen bulb and automatic timer
- Integrated laser beam for precise distance measurement
- ▶ 2-Speed manual exposure switch on 2.5 m curled cable
- Includes practical and sturdy aluminium case
- CE Certified & PTB approved and tested
- Dimensions (in cm): L 28 x W 16.3 x H 16
- Lightweight: 8.9 kg

701115



HIRAY 5 X-RAY MACHINE

The portable HiRay 5 combines ultimate power with extreme compactness. It produces a nominal output of 5 kW, while weighing only 14 kg! The greater mA range allows for shorter exposure times. This is ideal for imaging regions with rapid motion such as the small animal thorax.

Technical specifications

- ► High-voltage setting: 120 kV / 100 mAs (30 mA)
- ▶ High-current setting: 50 kV / 32 mAs (100 mA)
- kV setting: 40 to 120 kV (in 1 kV steps)
- ▶ 26 mAs settings, 0.3 100 mAs, with illuminated LCD display
- Generator: High frequency inverter, 5 kW nominal output
- X-Ray Tube Focal spot: 1.8 mm
- Collimator with bright halogen bulb and automatic timer
- Integrated laser beam for easy distance measurement
- 2-speed manual exposure switch on 2.5 m curled cable
- Dimensions (cm): L 35 x W 20.3 x H 16 cm
- Weight: 14 kg
- ▶ Foot switch optional

701105



TABLE STAND FOR MOBILE X-RAY DEVICES

- "Floating" tabletop considerably simplifies patient positioning (4-way tabletop with 60 cm longitudinal and 20 cm lateral displacement, mechanically braked, L 150 x W 72 x H 78 cm including 4 fixing clamps)
- Self-adjusting bucky tray, 45 x 45 cm
- ▶ Dimensions (in cm): W 95 x L 150 x H 200



EQUINE X-RAY SYSTEM

For digital X-ray imaging in the mobile equine practice

The DigiVet DR Equine wireless panel set impresses with its simple operation, fast, wireless transmission speed and razor-sharp, high-resolution X-ray images.

- Active area of 25 x 30 cm (10 x 12")
- ▶ Highest resolution through CsI PLUS scintillator
- ► Stable and reliable automatic exposure detection (Anytime™)
- DQE technology for reduced dose
- X-ray voltage range: 40 150 kVp
- ▶ Battery: Lithium Ion, 8 hours standby or approx. 750 shots
- Dimensions (in cm): L 28.7 x W 35 x H 1.5
- Weight: 2.1 kg

Set consists of:

- DigiVet DR Equine WiFi panel
- VXvue™ image editing software
- ▶ High-end notebook with 15.6" full HD monitor
- ▶ PODOBLOCK® protection box
- Outdoor case for all components

708500

Online Training S708482

HiRay LITE portable X-Ray Machine

The HiRay LITE is a compact, battery-operated X-Ray device for mobile equine practices. With a weight of only 7.2 kg and a power output of up to 90 kV at 20 mA, it is well ahead of most competitors in its weight class.

- ▶ High RF power: 50 90 kV / max 20 mA, 1.6 kW
- High resolution imaging at low dose
- Focus: 1.2 x 1.2 mm
- ▶ Bright LED collimator light
- Integrated dual laser beam for easy distance measurement
- 2-stage manual release switch
- Lithium-Polymer Battery for up to 560 shots
- Less than 3 hours charging time
- Compact dimensions (in cm): L 31 x B 20 x H 17
- Light weight: only 7.2 kg
- Incl. transport case and charger





708500



701111



DIGIVET DR FLAT PANEL DETECTOR

The DigiVet DR is the latest easy DR solution that requires no change to your conventional X-Ray System. A large coverage area of 17×17 " (43 cm x 43 cm) enables a quick and effortless integration. This solution not only enables ultra-fast acquisition times of only 1.2 seconds but also features an amazing 9.4 Megapixel resolution. The DigiVet DR comes with VXvue, a user-friendly and versatile imaging software that is compliant with the DICOM 3.0 standard.

Scintillator: CsI (TI)

Panel technology: a-Si TFT with PIN diode

Pixel Pitch: 140 x 140 um

Resolution: 9.4 Megapixel (3,072 x 3,072 pixels)

Image size: 17" x 17" (43 cm x 43 cm)

A/D Conversion: 14 bitGrayscale: 16,384 grayscalesInterface: Gigabet Ethernet

Automatic Exposure Detection (AED) imaging technology

Dimensions (in cm): W 47 x L 47 x H 1.5

Weight: approx. 4.5 kg

708482

Workstation PC with 23.8" EIZO Monitor

English version incl. QWERTY-Keyboard 708217



708482

X-RAY APRON

- Lightweight and comfortable protection
- Water repellent and stain resistant nylon
- Velcro side straps
- 0.5 mm lead protection
- Colour: Dark blue

Medium

Shoulder: 41 cm width

▶ Hips: 60 cm width

▶ Height: 170 - 180 cm

704935

Large

Shoulder: 43 cm widthHips: 60 cm width

Height: 180 cm



704935

PROFIDENT PLUS DENTAL UNIT

The Profident PLUS is the result of many years of experience. It offers a wide range of applications for both, general practitioners and dental specialists. It is characterised by its modern design combined with high mobility. The stand is adjustable in height and all handpieces are arranged within easy reach of the operator. When a handpiece is removed from its holder, the Profident PLUS automatically supplies air and water to it.

High Speed Turbine (300,000 rpm)

- Push button burr change for friction grip burrs (FG)
- Variable burr water cooling facility
- With LED Illumination

Low Speed Air Motor (20,000 rpm) with Straight Handpiece

- For long straight shank burrs (HP), diamond discs and prophy heads
- High torque, forward and reverse rotation
- Optional water cooling facility (item no. 174022)

Built-in Dental Scaler

- ▶ Handpiece with LED Illumination
- Includes 5 scaling tips

3 in 1 Syringe

- Precise finger tip control of water, mist or dry air
- Water connection to connect external dental scaler (e.g. item no. 176870)

Dental Unit Air Compressor

- Powerful (340 watt)
- Quiet, 40 dB (A) (as a fan)
- Low maintenance

Micro Oil Filter

- Reduces the oil content in the compressed air down to 0.01 %
- Helps extend the life of handpieces

174001

Starter Kit

Complete set, consisting of:

- ▶ 1x Bur Kit for Contra Angle
- 1x Prophy Cups
- 1x Polishing Paste
- 1x Cleaning Brush

EXP174004



Scaler with LED lighting



Dental Handpiece High Speed (FG)



174001

MOBILE DENTAL CART

Turns every table into a dental treatment table!

With the Mobile Dental Cart from EICKEMEYER®, every operating or treatment table becomes an ergonomic dental treatment table.

- ▶ For attachment to operating and treatment tables
- Ergonomic, curved dental tray with grate
- Easily removable tray for cleaning
- Made of stainless steel
- With 4 lockable smooth running castors
- ▶ Height adjustable from 75 100 cm
- Tray size (in cm): W 50 x D 35





COCOON SPRAY DENTAL UNIT

Its compact and robust design makes this unit ideal for the majority of all routine veterinary dental procedures.

The dental scaler is a state-of-the-art Piezo system, it offers superb efficiency to remove even the most stubborn tartar with ease.

Water cooled drilling means you can run faster speeds without the risk of thermal tissue damage.

- Infinitely variable power control
- Self-regulated frequency: 27 to 32 kHz
- Driver type: piezoelectric
- Oscillation form: linear
- For straight and angled dental handpieces
- With water cooling
- ▶ High torque to ensure maximum precision
- Self-regulated speed up to 40,000 rpm
- Handpieces not included

176860

Low Speed Dental Handpiece - Contra angle (RA)

- With water cooling
- Compatible with Cocoon Dental Unit (176860)
- For polishing please order: Polishing Cups RA (176787)

174037



EICKSONIC II DENTAL SCALER

- With illuminated LED handpiece (removable)
- High performance (28 kHz)
- With modern Piezo technologies
- Suction power control
- ▶ Water connection hose with connector for 3/8" water line
- Very easy to use
- Compact, durable design
- Best price-performance ratio
- Includes: 5 tips (2x flat/ wide; 2x flat/ pointed; 1x punctiform)



EICKEMEYER® DENTAL POLISHER / DRILL UNIT II

ORDER ONLINE www.eickemeyer.com

- Compact dental micro motor for straight and angled dental handpieces
- Powerful and durable, quiet operation
- Infinitely variable speed up to 35,000 rpm
- Forward and reverse movement
- Choice of hand or foot switch control
- Handpieces not included

174028



DENTAL STARTER KIT FOR SMALL ANIMALS

The most important equipment for Dental Unit Polisher / Drill Unit II (174028) and Dental Unit EickSonic Plus LED (176880) together in one set.

Complete set, consisting of:

- ▶ Low Speed Dental Handpiece Contra Angle (RA) (174026)
- Low Speed Dental Handpiece Straight (HP) (174027)
- Dental Burr Set (RA) (174105)
- Disposable polishing heads (HP) (174066)
- ▶ EICKEMEYER® Prophy Paste (176782)
- Dental Handpiece Lubricant Oil (174010)
- Dental Handpiece Air Motor / Turbine Adapter (17401150)
- ▶ Dental Handpiece Low Speed E-Type Adapter (17401155)

174002





174105

175782

DENTANOMIC DENTAL INSTRUMENTS

Easier, safer & faster dental surgery

DENTANOMIC instruments are premium surgical tools with superb tempered and hardened modular blades making dental extraction in dogs and cats easier, safer and faster.

- Blades are made from top quality, hardened stainless steel and cut much better than conventional instruments
- Unique modular design allows easy blade replacement when blades are blunt or damaged
- Ergonomic handle design promotes safe palm grip, only minimal finger force is required to stabilise the cutting edge
- Standard handles are perfectly balanced and sized for surgeons with hands up to 7.5 glove size



DENTANOMIC Dental Instruments	
Item no.	Description
175000	Ergonomic Handle, blue
175001	Ergonomic Handle, green
175002	Ergonomic Handle, pink
175008	Ergonomic Handle, orange
175009	Ergonomic Handle, purple
175005	Ergonomic Handle, blue, longer shaft
175006	Ergonomic Handle, green, longer shaft
175007	Ergonomic Handle, pink, longer shaft
175003	Ergonomic Handle, orange, longer shaft
175004	Ergonomic Handle, purple, longer shaft
175100	Elevator Tip, 2 mm
175101	Elevator Tip, 3 mm
175102	Elevator Tip, 4 mm
175103	Elevator Tip, 5 mm
175104	Luxation Tip, 1.5 mm
175105	Luxation Tip, 2 mm
175106	Luxation Tip, 3 mm
175107	Luxation Tip, 4 mm
175108	Luxation Tip, 5 mm
175200	Storage and Sterilisation Rack for 9 blades
175202	Sterilisation Container for 5 handles
175201	Honing Pack, extra-fine

DENTANOMIC PERIOSTEAL ELEVATOR SET

DENTANOMIC dental instruments are an innovation in the field of dental surgery and significantly simplify tooth extraction. With the interchangeable blades, the instruments can be quickly and easily adapted to the respective surgical situation. In addition, once the blade wears out, the entire instrument no longer needs to be thrown away. The blades are extra hardened and have a special, very sharp grind.

Set consisting of 4 periosteal elevator blades (2, 4, 6 and 8 mm) and sterilization container



DENTANOMIC DENTAL INSTRUMENTS

DENTANOMIC Set - Bronze

Complete Set, consisting of 2 handles (purple and gold) + elevator blade set (all sizes) + luxator blade set (all sizes) + sterilisation and storage rack

Up to glove size 7.5 From glove size 8.0

175050 175051

DENTANOMIC Set - Silver

Complete Set, consisting of 3 handles (green, blue and pink) + elevator blade set (all sizes) + luxator blade set (all sizes) + sterilisation and storage rack + honing set

Up to glove size 7.5 From glove size 8.0

175060 175061

DENTANOMIC Set - Gold

Complete Set, consisting of 5 handles (purple, gold, green, blue and pink) + elevator blade set (all sizes) + luxator blade set (all sizes) + sterilisation and storage rack + honing set + sterilisation container

Up to glove size 7.5 From glove size 8.0

175070 175071

DENTANOMIC Set - Gold Mixed

Complete Set, consisting of 3 standard handles (green, blue and pink – up to glove size 7.5) and 2 long handles (purple and orange – from glove size 8.0) + elevator blade set (all sizes) + luxator blade set (all sizes) + sterilisation and storage rack + honing set + sterilisation container




EICKWARM H₂O HEATING SYSTEM FOR SMALL ANIMALS

Prevent Intraoperative Hypothermia with the EickWarm H2O Patient Heating System

Intraoperative hypothermia is one of the most common anaesthesia complications in small animal medicine. Effective perioperative heat management is crucial for reducing this risk and ensuring patient safety.

The EickWarm H2O Patient Heating System provides safe and effective heat therapy with precise temperature control for small animals. By circulating demineralized water through a leak-proof heating mat, the system delivers conductive temperature therapy to prevent heat loss and hypothermia during surgery and recovery. The compact device operates quietly, ensuring a calm environment for both patients and staff. The efficient heating element quickly brings the warming mat to the desired temperature, providing immediate warmth. Additionally, by adding pre-cooled water, the device can be used to cool patients, offering versatile temperature management.

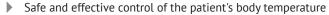

- Continuous temperature adjustment from 10 °C to 45 °C
- For use with demineralised water
- Digital temperature and water level display
- ▶ Three different timer modes
- Locking valves prevent the device from leaking when changing mats
- Includes two warming mats (L 55 x W 55 and L 55 x W 100 cm)
- Compact dimensions (in cm): L 19 x W 19 x H 17
- Weight: 1.2 kg

648020

Small Replacement Mat for EickWarm H₂O

Dimensions (in cm) L 55 x W 55

E64802001

Large Replacement Mat for EickWarm H₂O


Dimensions (in cm) L 55 x W 100

E64802002





E64802001

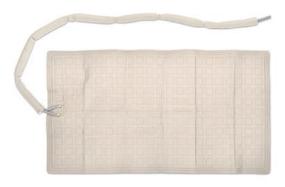

E64802002


REUSABLE WARMING MATS

These reusable Warming Mats are placed under the patient and – thanks to their unique 3D textile structure – ensure optimal heat supply. The mats fit the body's contour perfectly, making temperature control possible in any position and without interruption at the contact points.

- Warm air is distributed evenly under the patient
- Pressure relieving positioning
- ▶ Free operating field
- Reusable and washable at 60 °C
- Compatible with most air warming devices
- Made in Germany

Small, 50 x 60 cm Large, 115 x 60 cm

648110 648120

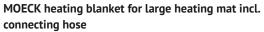

- Especially suitable for dental interventions
- In combination with the large heating mat, ensures optimal convective heat supply from above and below
- Variable round neck suitable for different head and body sizes

648130





648110



648110

EICKWARM AIR WARMING SYSTEM FOR SMALL ANIMALS

Intraoperative hypothermia is one of the most common anaesthesia complications in small animal medicine. Effective perioperative heat management is crucial for reducing this risk and ensuring patient safety. The new EickWarm Air Warming System quietly delivers forced air through optionally available warming mats that are specifically designed to evenly deliver diffused warm air to the patient.

- Safe and effective management of patient body temperature
- Optimal performance through uniform heat distribution
- Extremely powerful and quiet (≤55 dB)
- Only 60 seconds to warm up
- Temperature set points: Ambient air, 32 °C, 35 °C, 38 °C and 43 °C
- ▶ Three selectable fan levels (air flow)
- Reliable measurement of the hose outlet temperature using sensors
- Multiple alarm options ensure maximum patient safety
- Pole clamp for secure fixation
- ▶ Long heat hose (2 m) with flexible fastening clamp
- Compact dimensions (in cm): L 30 x W 37 x H 34
- Weight (including heat hose): 6.5 kg
- Warming mats available separately



EICKLOXX MICRO OSTEOSYNTHESIS SYSTEM

The EickLoxx Micro osteosynthesis system was specially developed for the treatment of fractures in very small dogs and cats.

The angle-stable screws have a diameter of 1.2 mm. In addition, various 1.0 mm thick plates are available in different configurations.

This opens up the following indication areas, especially for dogs and cats weighing less than 3 kg as well as puppies: radius/ulna fractures, metatarsal fractures, jaw, pelvic and scapula fractures, in addition to arthrodesis and much more.

The EickLoxx Micro osteosynthesis system joins the existing EickLoxx osteosynthesis system for dogs and cats and supplements this in the increasingly important lowest weight category.

197400



EICKLOXX SMALL CONTOUR OSTEOSYNTHESIS SYSTEM


The EickLoxx Small Contour plate set includes various anatomically preformed osteosynthesis plates for the treatment of fractures of the distal radius and femur, ilium, acetabulum and for arthrodesis of the carpus and tarsus. It complements the cut-to-size plates of the EickLoxx Small osteosynthesis system and is therefore also designed for the treatment of cats and dogs up to 15 kg. The instruments from the EickLoxx Small osteosynthesis system are used in a similar way.

The titanium plates of the EickLoxx Small Contour plate set are available in two plate thicknesses: The 2.2 mm plates (yellow, green and magenta) can be locked polyaxially with the 1.7 mm and 2.3 mm screws and the 3.0 mm plates (light blue) optionally with the 2.7 mm, 3.5 mm and 4.0 mm screws of the EickLoxx series.

Despite contouring, the panels can be bent as required using the cabinet iron and / or the panel bending pliers with rollers and shortened using the cutting pliers.

The sieve tray with the EickLoxx Small Contour plates fits into the container Art. No. 185555 on the two sieve trays of the EickLoxx Small system.



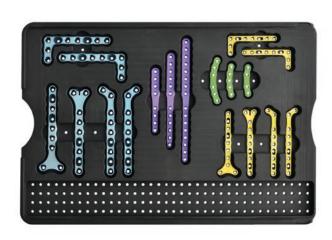

EICKLOXX SMALL OSTEOSYNTHESIS SYSTEM




The system distinguishes itself through the polyaxial placement of the custom-fit screws and combines the advantages of a locking plate system with polyaxial screw placement in $\pm\,15^\circ$ longitudinal and transverse pivoting.

Compared to conventional osteosynthesis systems, the use of locking systems (internal fixation) such as the EickLoxx Small facilitates faster healing of fractures and improved implant strength. It is for this reason that the removal of plates and screws is rarely indicated.

The biologically advantageous titanium plates are available in three sizes and can be cut to the desired length and bent in three planes. Specially designed tools are included to perform this. Titanium also reduces artefacts in post-operative imaging.

EickLoxx Small is a modular system. All implants and instruments are made in Germany using only the highest quality materials. **185500**



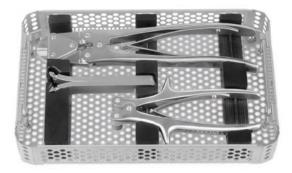




EICKLOXX LARGE OSTEOSYNTHESIS SYSTEM

The universal polyaxial osteosynthesis system for small animals The newly developed EickLoxx Large is a polyaxial locking system for dogs with a weight of approx. 15 to 45 kg.

Just like the EickLoxx Small, the system is characterised by a polyaxial placement of the precisely shaped screws and therefore combines the advantages of locking systems and the possibility of placing screws at an angle of up to ± 15°. Like the EickLoxx Small, the EickLoxx Large offers the advantage of rapid fracture healing with improved implant strength. It is for this reason that the removal of plates and screws is rarely indicated.

The biologically advantageous titanium S-shaped plates are available in eight different sizes. Thanks to special tools, the plates can be bent on all three planes. The long 12- and 14-hole plates as well as the long 34 to 40 mm screws (Ø 2.7 / 3.5 / 4.0) are not included in the basic set and can be added if required. Both plates can be integrated into the implant tray. There is space in the implant tray for the long plates and capacity for the module for the long screws in the instrument tray.

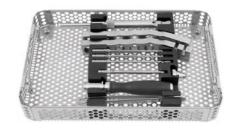





ZLIG - NEW TECHNIOUE FOR CRUCIATE LIGAMENT REPLACEMENT (CRCL)

Zlig - the history

The tear of the anterior cruciate ligament is still one of the most common orthopaedic diseases in dogs. The path of the many surgical methods developed for this vary from intracapsular and extracapsular techniques to modern corrective osteotomies that alter the geometry of the affected knee joint. With the development of new materials in medical technology, it is now possible to replace the cranial cruciate ligament in an anatomically correct manner, instead of changing the forces acting on the joint. After a long period of preparatory work by the French Dr. Jacques-Phillipe Laboureau, a suitable synthetic ligament is now available for the intra-articular cruciate ligament replacement in small animals. Together with the instrumentation developed by EICKEMEYER®, this new technique for cruciate ligament replacement can now be performed.

Zlig - the implant

The Zlig consists of ultra-high-molecular polyethylene with the special feature that the woven structure of the implant is interrupted intra-articularly by "free fibres". Free parallel fibres reduce fatigue and also encourage the ingrowth of fibroblasts and collagen. Each implant is delivered in a sterile packed sleeve, making it easier to handle the implant before it is inserted into the joint. A selection of sizes with different resistances and fibre lengths are available to fit different patient sizes.

Zlig - the technology

In this technique, an artificial ligament is used as a total replacement for the anterior cruciate ligament using tunnel-tunnel-technique. The ligament is fixed in the tibia and Os femoris using specially developed cannulated interference screws in drill channels. The screws are guided parallel to the ligament using a guiding wire to avoid deviations. The technology does not cause irreversible damage and the results are reproducible thanks to a quick learning curve. Another great advantage of the technique is the fact that patients can put weight on the hind leg immediately after the operation without any risks.

16 Fibres / 10 mm

▶ 5-8 kg

2,000 N

191501

24 Fibres / 15 mm

▶ 8-12 kg

3,000 N

191502

32 Fibres / 17 mm

▶ 12-25 kg

▶ 4,000 N

191503

48 Fibres / 19 mm

> 25 - 45 kg

6,000 N

191504

48 Fibres / 22 mm

> 25-45 kg

▶ 6.000 N


191505

48 Fibres / 25 mm

> 25-45 kg

▶ 6,000 N

191506

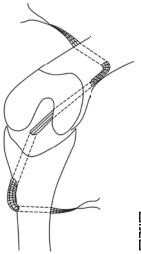



QR code for the application video of the Zlig

ZLIG 2 – NEW TECHNIQUE FOR CRUCIATE LIGAMENT REPLACEMENT (CRCL)

Zlig 2 - the instrumentation

A small and inexpensive set of instruments is required to perform this innovative cruciate ligament surgery method. The threads of the specially developed titanium interference screws are round so they do not cause any damage to the fibres of the Zlig.





ORTHOVET PRO BATTERY POWERED DRILL SYSTEM



The ideal tool for small animal osteosynthesis

The OrthoVet PRO is a modular, battery driven drill system designed for small animal osteosynthesis. The modular system consists of a handpiece and four attachments, which can be changed intraoperatively in just a few seconds with a keyless quick coupling.

The handpiece, the quick release AO Chuck and the Jacobs Chuck are cannulated. Thus, cannulated screws or drills can also be tracked over long guide wires in stages, without the guide wire having to be shortened. The oscillating saw has a fivefold adjustable, variable fixation of the saw blade and allows the surgeon an optimised view of the surgical field. The AO quick coupling completes the system and makes the OrthoVet PRO a real all-rounder.

The handpiece and attachments can be sterilised / autoclaved. The Li-Ion Battery is not autoclavable and is inserted into the handpiece by means of special holder. Continuous operation is possible in clockwise and counterclockwise motion, as well as in oscillation with the saw attachment.

- 4 attachments: Jacobs Drill Chuck, Quick Release AO Chuck, Quick Release Chuck for KIRSCHNER wires, Oscillating Saw with Saw Blades
- ▶ Fully Cannulated Handpiece
- Autoclavable
- Includes 2 Li-Ion Batteries with Charger and Transport Case

Technical specifications

Infinitely variable speed: 0 – 1,200 rpmOscillation mode: 0 – 16,000 ocs/min

▶ Battery type: Li-Ion

Battery charging time: 180 minWeight (includes Battery): 933 g





The powerful and intuitively operated motor system for small animal osteosynthesis

The DRAGONFLY is a modular cordless drill system with 4 attachments for screwing, milling, drilling and cutting as required. Thanks to the quick coupling, the attachments can be changed intraoperatively. The handpiece itself impresses with its narrow handle design and an intuitive one-handed mechanism.

The AO attachment is cannulated, same as is the keyless chuck for KIRSCHNER drilled wires / STEINMANN nails and the Jacobs drill chuck. Both are equipped with an adjustable clamping range of up to 4.0 mm. These 3 attachments can be operated stepless in left / right rotation and also in oscillation mode at speeds from 0-1,350 rpm.

The DRAGONFLY is completed by attachment no. 4: the oscillating saw. It is equipped with a power of 0–14,000 osc/min and allows the saw blade to be used in 4 different positions. Handpiece and attachments can be autoclaved. The core, however, is the removable, non-autoclavable drive unit: This includes the control unit, the battery, the electronics and the motor control electronics and, together with the Li-Ion technology, guarantees a long service life and service life. Thus, the DRAGONFLY offers you a performance that you previously only knew from compressed air or cable-electric systems. All components are designed and manufactured in Germany.

Technical specifications

- Infinitely variable speed: 0 1,350 rpm
- Oscillation cutting mode: 0 14,000 osc/min
- Battery type: Li-lon
- Scope of delivery: DRAGONFLY handpiece incl. 4 attachments (AO attachment, keyless chuck for KIRSCHNER bone wires / STEINMANN nails, Jacobs drill chuck, oscillating saw), 4 saw blades with STRYKER attachment, drive unit, sterile funnel, charger and carrying case











ICARE® TONOVET PLUS TONOMETER

The Icare® TONOVET Plus is easy to use. It features an advanced user interface with a clear, large display. Red/green positioning lights help to find the correct angle and aim at the middle of the cornea. The automatic measuring sequence lets you choose between taking single readings or taking all 6 measurements with one press of the button.

- Quick, painless and accurate procedure
- Easy to use and requires no calibration
- Assisting positioning lights
- Large and easy to read colour display
- ▶ Dimensions (in mm): L 215 x W 24-29 x D 35-95
- Weight: 230 g (incl. batteries)
- Includes: Transport Case, 100 Probes, 4x AA batteries

1736007



QR code for the application video of the TONOVET Plus

Stand for TONOVET and TONOVET Plus

Table and wall mountable

1736008





Application



ICARE® TONOVET PRO TONOMETER

The iCare TONOVET Pro tonometer sets new standards in intraocular pressure (IOP) measurement and is aimed at vets who have the highest demands on accuracy and ease of use. The innovative "Quick Measure" function enables faster and easier measurement of the IOP and significantly shortens the examination time, while the more responsive test tip and the optimised test tip control ensure improved handling and even more precise measurements. The 200-degree positioning freedom makes the TONOVET Pro tonometer an indispensable helper even in emergency situations. The ergonomic display is designed to allow the user to instanty read the measured values.

- Improved user-friendliness
- 200-degree positioning freedom (idela for equine use)
- More responsive test tip
- Better test tip control improves accuracy
- Display facing the user
- Self-calibrating
- Wireless printing via Bluetooth
- ▶ Includes: Transport Case, 100 TONOVET Probes, 4x AA batteries

173625

Probes for TONOVET

▶ 100 per pack



Application





173625

KOWA SL-19 PORTABLE SLIT LAMP

Now more ergonomic and even lighter!

With its rounded shape and curved handle design, the new KOWA SL-19 Slit Lamp is even more ergonomic. Furthermore, the weight of the hand-held unit has been significantly reduced.

- New LED backlight that can be used simultaneously with the slit light
- New LED blue light illumination for flourescence observation
- Light intensity can be adjusted with one hand
- Infinitely adjustable slit length (1.5 12 mm)
- Longer operating time: up to 6 hours at maximum intensity
- Available in 5 colors: white, black, blue, green and pink

Technical specifications

- Magnification: 10x and 16x
- Working distance: about 80 mm at 16x magnification, about 100 mm at 10x magnification
- ▶ Effective field of view: Ø 10 mm at 16x magnification, Ø 15 mm at 10x magnification
- ▶ Slit width: Ø 0.1, Ø0.2, Ø 0.8 mm
- Spot size: Ø 1, Ø 5, Ø 12 mm (slit length continuously adjustable from 1.5 − 12 mm)
- Light intensity: 200 20,000 lux (continuously adjustable)
- ▶ Power supply: 3x AA batteries (DC 4.8 6.4 V)
- Dimensions (in mm): W 197 x D 107 x H 238
- Weight (without batteries): only 620 g

173690

Carrying Case for SL-19 173695



173690

KOWA SL-19 PLUS PORTABLE SLIT LAMP

The latest development from KOWA, the hand-held slit lamp SL-19 plus, impresses with its well-known ergonomic shape and light weight with an integrated, high-resolution camera. With a simple hand movement, videos can be recorded during the examination, which are automatically transmitted to an iPad (not included) via WiFi. Using special application software, these videos can be analysed and further processed on the iPad. High-resolution still images can also be extracted from the videos, which can be individually selected, saved and forwarded.

Technical specifications

- Magnification: 10x
- Working distance: approx. 100 mm
- ▶ Effective field of view: Ø 15 mm at 10x magnification
- Camera field of view: approx. 27 to 15 mm
- ▶ Effective pixel count: approx. 2 megapixels
- A camera indicator shows the operating status of the camera
- Operating time: approx. 8.5 hours if the device is switched off and disconnected from the iPad after each use, approx. 3.5 hours at full power without disconnecting from the iPad
- Dimensions (in mm): D 197 x W 107 x H 255
- Weight (without batteries): 730 g

173696

Carrying Case for SL-19 plus 173697



173696

BAUSCH & LOMB PUREVISION® 2 BANDAGE CONTACT LENSES

ORDER ONLINE www.eickemeyer.com

Contact bandage lenses act as a protective shield on the cornea. They permit corneal repair and improve patient comfort by preventing contact between damaged cornea and eyelids.

- ▶ High definition optics offers a clear and sharp view
- Excellent stabilisation
- ▶ Comfortable contact lens wear for up to 30 days and nights
- Ensures an optimal supply of oxygen to the cornea and makes eyes look clear
- Packaging: 6 lenses per box

173540



ALGERBRUSH II

The AlgerBrush II is a battery powered, handheld, low torque motorised burr. It has been designed to remove foreign bodies from the cornea or sclera of patients.

The burr in this rotary instrument is used as a "brush" rather than a drill and the foreign bodies are thus "brushed" from the eye leaving a smooth surface, which in turn heals much faster than with other instruments that could be used for this purpose.

According to the paper "Veterinary Ophthalmology" ([2013] 16, 2, [83-88]), it has been successfully used to treat spontaneous corneal epithelial defects (SCCEDs) in dogs in combination with a bandage contact lens.

Includes

- > 3.5 mm round, medium grit diamond burr
- Chuck

170010

AlgerBrush II Corneal Burr

- Ø 3.5 mm
- Chuck attachment included

170011

AlgerBrush II Corneal Burr

Ø 2.5 mm

170013

AlgerBrush II Corneal Burr

- Ø 1 mm
- Pack of 5 pieces

170012

AlgerBrush II Diamond Grinding Wheel (4 mm)

Fine



170010









1/0014

HEINE OMEGA 600 INDIRECT OPHTHALMOSCOPE

Discover more. With less weight.

Lighter. Brighter. Better. HEINE kept the best and perfected the rest. Highly comfortable, the HEINE OMEGA 600 delivers a new class of technological innovation, along with the new and unique visionBOOST for significantly better diagnostic even in dense cataracts.

- Most precise diagnostic, up to the far retinal periphery
- ▶ Up to 20 % better view in media opacities with visionBOOST
- Lightest High-End Indirect Ophthalmoscope worldwide
- Long-lasting investment due to its robust design
- Easy disinfection thanks to smoothSURFACE

Most precise diagnostic, up to the far retinal periphery

- Significantly better diagnostic due to unique visionBOOST
- Sharpest fundus image due to superior glass optics
- Brightest fundus image with true colour rendering thanks to our unique LEDHQ – now in even more natural colours
- ▶ Best 3D imaging due to our unique Stereoscopic Adjustment
- Technology in dilated and undilated pupils

Includes

OMEGA 600 Indirect Ophthalmoscope (with battery and +2 D eyepiece), E4-USBC power supply with cable, eyepiece 0 D

173497

CB1 Charging Battery 17349702



173497

HEINE DIAGNOSTIC SET

Absolute precision is the basic requirement to ensure the creation of a high-quality product. At HEINE this principle applies to the entire production process.

Complete set, consisting of:

NT 300 Double Charger (302550)

- Automatic charging process and overcharge protection
- Includes two NT 3.5 V Battery Handles including batteries (2x 302535)
- With short-circuit protection and a continuous ceramic brightness control

BETA 200 LED Ophthalmoscope (300137)

The HEINE BETA 200 LED Ophthalmoscope offers a bright, clear view with durable LEDHQ technology and precise dimming, ideal for reliable mobile diagnostics.

G100 LED Slit Illumination Head (300036)

The HEINE G100 VET LED features 1.6x magnification and precise, dimmable LED HQ lighting in a robust all-metal housing, ensuring optimal visibility and reliable diagnostics.

Replacement lamp can be ordered through item no. 300038.

Otoscope Specula

- High quality chrome plated metal construction, which is easy to clean, disinfect and sterilise
- Metal Speculum, 4x 57 mm (300357)
- Metal Speculum, 6x 65 mm (300365)
- Metal Speculum, 7x 90 mm (300390)

EXP302550



HEINE LED LOUPELIGHT 2 WITH 2.5X BINOCULAR LOUPES



LED in HQ – the new standard in LED lighting delivers homogeneity and bright light with natural colour rendering for a precise diagnosis. The combination of light and magnification provides optimal results.

HEINE LED LoupeLight 2

- Illumination: 55,000 lux
- Colour temperature: 4,000 Kelvin
- Adjustable viewing angle ensures a coaxial, shadow-free image
- Replaceable connecting cable, 1.6 m
- Including HEINE mPack mini and USB plug-in power supply
- Operating time ~ 9 hours, at full power ~ 4 hours

310717

HEINE HR binocular loupes with 2.5x magnification

- With i-View loupe mount on S-FRAME
- Including protective glasses, 2 sterilisable swivel levers, retaining strap, cleaning fluid, cleaning cloth and case
- Working distance: 340 mm
- Field of view: 90 mm
- Depth of view: 110 mm

310700

- Working distance: 420 mm
- Field of view: 130 mm
- Depth of view: 180 mm

310707



HEINE ML4 LED HEADLIGHT WITH MPACK UNPLUGGED & HEINE HR LOUPES

Perfect fit. Perfect view. Even with long operations.

The HEINE ML4 LED HeadLight is specially designed for healthcare professionals who need a perfect view during long surgical operations or treatments. Various individual adjustment options of the headband and the comfortable leather finish guarantee the highest level of comfort and a secure fit.

HEINE ML4 LED HeadLight with mPack UNPLUGGED

- Lightweight, high-performance battery integrated on the back of the headband
- Offers 100% freedom of movement
- Incl. Plug-in Transformer UNPLUGGED

311138

HEINE HR Loupes

- With i-View loupe holder
- ▶ 3.5x magnification
- ▶ 420 mm working distance

311123

- 4x magnification
- > 340 mm working distance



LED DIAGNOSTIC LIGHT

This high quality diagnostic light brings light into the dark for a wide range of applications – whether used with or without the spatula.

- Ideal for intubating small animals
- Convenient for examining eyes, ears, nose and throat
- ▶ Helpful for dental treatments, especially on rodents
- White LED with a lifespan of more than 100,000 hours
- Incl. 1 angled spatula and 2 AAA batteries

214200



EICKVIEW COMPACT ENDOSCOPE SET

The specially developed EickView Compact Endoscope provides the perfect solution for rhinoscopy, tracheoscopy and vaginoscopy in small animals due to its small diameter (3.8 mm) and integrated working channel.

- Highly mobile
- Homogeneous and minimally invasive distal end
- ▶ Integrated working channel Ø 1.3 mm
- Outer diameter Ø 3.8 mm
- Working length 14 cm
- Viewing angle: 30°
- Opening angle 90° (wide-angle)
- Foreign body forceps Ø 1 mm, L 25 cm
- Attachable LED Light source
- Long battery life and spare battery with an operating time of 4 hours each
- 3 Adjustable light levels
- Stopcocks: two integrated stopcocks allow steady insufflation and irrigation
- ▶ Connection: STORZ, WOLF or ACMI connectors
- Eyepiece: DIN ocular allows the connection of all standard C-mount cameras
- ▶ Set includes: EickView Compact Endoscope incl. fitted otoscopy protection aid, portable LED light source incl. charger, foreign body forceps, spare battery, 2 cleaning brushes, padded carrying case

306122

Endoscopy Biopsy Forceps

- For taking biopsies
- Length: 250 mm
- ▶ Ø 1.0 mm



EICKVIEW SD ENDOSCOPY SYSTEM



Unlock the future of veterinary care with the EickView SD Endoscopy System!

This state-of-the-art system features a high-quality endoscopy monitor with a large, intuitive touchscreen, seamlessly connecting to three versatile "semi-disposable" endoscopes. This system allows the veterinary practitioner to experience the benefits of flexible endoscopy (Gastroscopy, Bronchoscopy, Rhinoscopy, etc.) without the hefty investment of traditional video endoscopes. Each endoscope is equipped with a high-resolution CMOS image sensor and an LED light source, ensuring crystal-clear visuals.

EickView SD Endoscopy Monitor

- ▶ 15.6" touchscreen with 1920 x 1080 Full HD resolution
- ► Functions: White Balance, Zoom, Freeze, photo (JPG) and video recording (MP4)
- ▶ 64 GB storage capacity
- Mains and battery operation, battery life up to 5 hours
- Outputs: HDMI, DVI, SDI and USB
- Dimensions (in cm): W 39 x H 25 x D 4
- Weight: 2.7 kg

307000

EickView 150SD

- Working length: 150 cm
- Outside Ø: 6.0 mm
- Working channel Ø: 2.9 mm
- Direction of view: 0°
- Viewing angle: 120°
- Bending up / down: 180°
- Bending right / left: 120°

307010

EickView 60SD

- Working length: 60 cm
- Ø outside: 3.2 mm
- Ø working channel: 1.2 mm
- Direction of view: 0°
- Viewing angle: 120°
- Bending up / down: 210°

307020



EickView 60SD-XS

- Working length: 60 cm
- Ø outside: 2.7 mm
- Ø working channel: 1.2 mm
- Direction of view: 0°
- Viewing angle: 120°
- Bending up / down: 210°



EICKVIEW PORTABLE LED VIDEO ENDOSCOPES

Due to their extremely compact dimensions, the EickView LED Video Endoscopes are the perfect solution for broncho-, gastro-, rhino-, tracheo, uretero- and vaginoscopy in dogs and cats and for the examination of the upper respiratory tract in horses. The devices work independently, as the powerful LED light source and the high resolution camera of the stand-alone systems are integrated in the distal end of the endoscopes. In addition, the folding 3.5" monitor allows both direct video and image playback as well as storage via SD card. Furthermore the EickView LED Video Endoscopes can be connected to a separate monitor.

- ▶ Folding 3.5" monitor for video and picture playback
- ▶ Recording and archiving on integrated SD card
- File formats: JPG, AVI
- Integrated powerful LED light source
- Rechargeable Battery with up to 3 hours of operation
- USB port and video out
- ▶ Depth of field: 5 60 mm
- Direction of view: 0°
- Field of view: 110°

EickView LED Video Endoscopes					
Item No.	Description	Working length	Outer Ø	Working channel	Deflection
306009	EickView 60S	60 cm	3.3 mm	1.2 mm	up 180° down 130°
306011	EickView 60	60 cm	4.2 mm	1.3 mm	up 180° down 130°
306010	EickView 70	70 cm	6.0 mm	2.0 mm	up 180° down 130°
306166	EickView 100E	100 cm	8.0 mm	2.8 mm	up 150° down 130°
306164	EickView 150E	150 cm	8.0 mm	2.8 mm	up 130° down 130°

Accessories				
Item No.	Description	Diameter	Length	
306211	Alligator Forceps, EickView 60/60S	1.0 mm	120 cm	
306311	Biopsy Forceps, EickView 60/60S	1.0 mm	120 cm	
306214	Alligator Forceps, EickView 70	1.8 mm	120 cm	
306314	Biopsy Forceps, EickView 70	1.8 mm	120 cm	
306216	Alligator Forceps, EickView 100E/150E	2.4 mm	230 cm	
306316	Biopsy Forceps, EickView 150E	2.4 mm	230 cm	



EICKVIEW 5-IN-1 ENDOSCOPY SYSTEM

EickView 5-in-1 Endoscopy System

The EickView 5-in-1 endoscopy system is a true all-rounder for veterinary medicine. The easy-to-use, mobile device offers the option to work with both, flexible and rigid endoscopes. This means that all types of endoscopy can be covered with just one device – from diagnosis to surgical procedures.

The EickView 5-in-1 combines a large, 15.6" Full HD display, a powerful LED light source, a connection for flexible endoscopes from the EickView HD series with air and water pump as well as digital data storage via SD Card in one system.

306200

EickView HD Video Endoscopes

With improved video chips, the EickView HD Video Endoscopes deliver even more brilliant images. The main areas of application are bronchoscopy and gastroscopy. Due to the large number of available scope sizes, they are suitable for all applications in small animals and equine.

Technical specifications


- Viewing direction: 0°
- Distal end contains 1 CMOS Camera Chip, 2 Light Guides, Working Channel, Lens Rinsing Nozzle
- Scope of delivery: EickView Video Endoscope, Pressure Tester, 2 Cleaning Brushes, Cleaning Set, Carrying Case

EickView 5-in-1 Camera Kit

The EickView 5-in-1 camera with zoom adapter and light guide cable completes the comprehensive package. It allows all standard rigid endoscopes to be connected to the main system.



306201

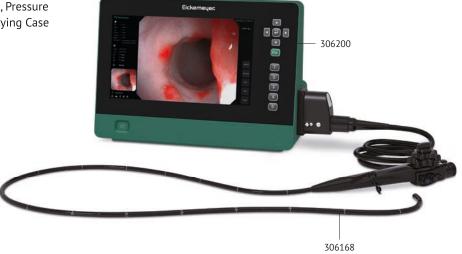

EickView HD Video Endoscopes							
Art. Nr.	Description	Working length	Outer Ø	Working channel	Depth of field	Field of view	Deflection
306177	EickView 60 HD Endoscope	60 cm	3.5 mm	1.2 mm	3 – 50 mm	120°	Angle above: 160° Angle below: 130° Angle right / left: 0°
306178	EickView 100 HD Endoscope	100 cm	6.0 mm	2.0 mm	3 – 50 mm	120°	Angle above: 210° Angle below: 90° Angle right / left: 100°
306168	EickView 150 HD Endoscope	150 cm	8.5 mm	2.0 mm	5-100 mm	≥ 140°	Angle above: 210° Angle below: 90° Angle right / left: 100°
306169	EickView 150L HD Endoscope	150 cm	9.2 mm	2.8 mm	5-100 mm	≥ 140°	Angle above: 210° Angle below: 90° Angle right / left: 100°
306170	EickView 300 HD Endoscope	300 cm	13.0 mm	3.2 mm	5-100 mm	≥ 140°	Angle above: 180° Angle below: 180° Angle right / left: 160°

ENBIO S B-CLASS AUTOCLAVE

Simple - reliable - fast

The Enbio S is the fastest B-class autoclave on the market. Thanks to its compact size, it fits into every practice and the sterilisation chamber offers enough space for a standard surgical kit. The Enbio S is set up in no time and it is simple and intuitive to use – switch it on, select the mode, sterilise – done!

- Extremely short program times of only 7 to 15 minutes
- Revolutionary low operating costs thanks to Magic Filter
- Very quiet operating noise of only 32 dBA
- Automatic documentation on USB stick
- ▶ Highest quality at a low price
- External dimensions (in mm): W 252 x H 162 x D 561
- Weight: 15 kg

500000



STERIHERO VET MELAG B CLASS AUTOCLAVE

Your new superpower in veterinary medicine.

No matter what big tasks the day throws at you, SteriHero Vet has your back – from the functions to the price! Powerful technologies, an intuitive smart touch display and the unique MELAG quality – made in Germany – give you the superpower to change everything!

- B Class autoclave for the safest sterilisation
- ▶ High loading capacities and short cycle times from 16 minutes
- Intuitive operating concept with large smart touch display
- Easy to fill, integrated water tank
- Holder for holding 5 trays (3 included)
- Chamber volume: 23 l, D 45 cm, Ø 25 cm
- Outside dimensions (in cm): W 47 x H 50 x D 64
- Weight: 49 kg



EICKSONIC ULTRASONIC CLEANERS

Universal ultrasonic cleaners with stainless steel tubs and stainless steel housing – for easy and safe instrument cleaning and excellent cleaning results.

- ▶ Timer: 0 30 min
- ▶ Heating: 20-80 °C
- Stainless steel basket protects cleaning material and vibrating tub
- Stainless steel lid reduces noise
- Operating frequency: approx. 40 kHz
- Mains voltage: 50 / 60 Hz
- Mains connection: 220 240 V
- Incl. stainless steel lid and screen basket

EickSonic 3.2

- Tub capacity: 3.2 liters
- Working capacity: 2.8 liters
- Ultrasonic power: 120 watts
- Heating capacity: 150 watts
- Without drain tap
- Dimensions:

Device – outside (in cm): L 26.5 x W 16.5 x H 22.5 Tub – inside (in cm): L 24 x W 13.5 x H 10

561032





EickSonic 6.5

- ▶ Tub capacity: 6.5 liters
- Working capacity: 5.8 liters
- ▶ Ultrasonic power: 180 watts
- ▶ Heating capacity: 300 watts
- With drain tap
- Dimensions:

Device – outside (in cm): L 32.5 x W 17.5 x H 27 Tub – inside (in cm): L 30 x W 15 x H 15

561065

EICKEMEYER® DOPPLER

Reliable and accurate non-invasive blood pressure monitor for small animals.

- Compact, hand-held design
- High sensibility for a clear signal
- Adjustable sound volume and headphone connection
- Focused 8 MHz Probe with wide beam for better flow detection in the artery
- Probe adapted for positioning on the animals tail or paw
- Rechargeable battery powered, operating time: ±6 hours
- Dimensions (in mm): L 175 x W 80 x H 25 mm
- Lightweight (300 g)



EICKEMEYER® MICROSCOPE

- A basic microscope with good optical quality at an affordable price
- Binocular head Siedentopf type, 30°
- ▶ Widefield eyepieces WF10x/18 mm
- Achromatic objectives 4x/0.1, 10x/0.25, 40x/0.65, 100x/1.25 oil
- Köhlersche illumination system, 12V /20 W halogen lamp
- Blue filter
- Abbe condenser 1.25 N.A.
- Mechanical stage 190 x 140 mm

710061



EICKEMEYER® LABORATORY CENTRIFUGE

- Safe and easy to operate
- Excellent value for money
- Holds 6 x 15 ml tubes
- Variable speed: 0 to 4500 RPM
- ▶ Built in timer, running time up to 60 min
- Safety latch stops the motor when the lid is open

718200



EICKEMEYER® REFRACTOMETER

- The refractometer is a valuable practice tool providing fast results
- Graduated scale
- Small sample sizes (1 drop)
- Lens with diopter adjustment




INFUSION STAND

- ▶ Height adjustable from 110 200 cm
- Stable plastic cross with 4 hooks
- Chrome finish
- ► Holder with dropper glass (plastic)
- 5 star castor base
- 2 brake castors

573106



SURGICAL DRAPES PVC

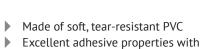

- Individually packaged sterile
- Pack of 25 pieces

spray adhesive

30 x 45 cm	90 x 120 cm	120 x 180 cm
920145	920212	921218
60 x 90 cm	120 x 120 cm	120 x 250 cm
920190	921212	921225



EICKWRAP FUN COHESIVE BANDAGE

The tried and tested EickWrap bandages are now also available with fun, colourful motifs. The bandages have a self-adhesive surface structure so that they stick to one another without any problems. EickWrap does not constrict and does not hinder blood circulation or desired joint movements. There is no sticking to skin, hair or clothing. Ideal for horse legs, hooves or dog and cat paws.

- Individually Packaged: Convenient and hygienic
- Length: 4.5 metres per roll / Width: 5 cm
- 6 Colourful Motifs
- ▶ Highest Quality Cohesive Bandage: Reliable and durable
- ▶ Conforms to Shape: Adapts to the area where it is applied
- Breathable: Allows air circulation over the skin and around the dressing
- Contains latex
- Box of 36 pieces





EICKTRON ELECTROSURGICAL UNIT

The EickTron is a versatile tool in veterinary medicine allowing it to be used for both monopolar and bipolar instruments. The unit's power adjustment via a foot or hand switch and its microprocessor control makes it very easy to use.

Includes

- ▶ Handpiece with finger switch and connecting cable
- ▶ Electrode set (10 pieces)
- Neutral electrode incl. cable
- Footswitch

Technical specifications

- ▶ Input: 220 240 V, 50 60 Hz
- Power consumption: 300 VA
- Frequency: 600 kHz
- Dimensions (in mm): L 288 x W 254 x H 104

323135



Bipolar Forceps

Angled

Length: 150 mm

323105

Bipolar Forceps Cable

323106

Bipolar Adapter

323139



CAT SCALE

- Contoured weighing surface
- Battery operated
- $\sim 0-20 \text{ kg}$
- Easy to clean
- ▶ Hold function
- Tare weighing function
- LCD display
- Button-on technology
- 5 g graduation (0.011 lb / 0.018 oz)
- 20 kg capacity (44 lb / 704 oz)
- Automatic switch-off function
- Toggle between kg / lb / oz
- ▶ Batteries included (2 x AA / 1.5 V)
- Dimensions (in cm): L 55 x B 31 x H 5.3



NOTES	
	_
	_

NUIES			

09/2024 | Photo titel: © Eric Isselée – stock.adobe.com