
veterinary technology for life
Eickemeyer®

Electrosurgery Unit EickTron 160 & Radiosurgery Unit EickTron 4.0

Item no. 323135 & 323400

The tips that are included with the handpiece are pictured on the right. They are all multi use. The handpiece should last for multiple uses with proper care. The flat rectangular tip that comes with the handpiece is the only tip that cannot be bought on its own – you get a new one with each new handpiece.

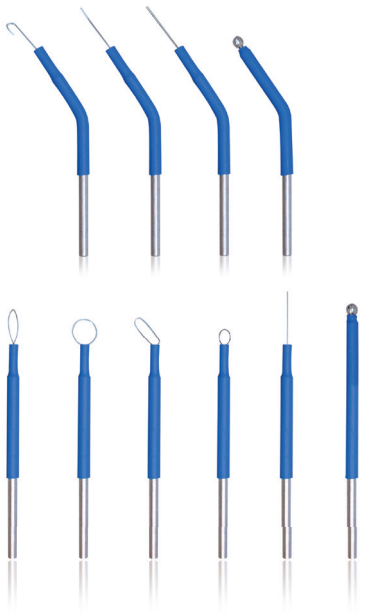

The handpiece as pictured on the right is steam autoclavable in a peel pack. The end tip also detaches, and we often recommend that this is kept in a small peel pack inside the other peel pack – that way when your RVT opens the handpiece if the tip hits the floor it is still sterile. You should be able to get multiple uses after autoclaving before you start to notice it not working as well.

Another trick that we recommend is always wrapping the metal grounding plate in a damp towel to help with connectivity to the patient.

For the settings, we typically suggest starting as low as 20 and then work your way up to a level that works. The lower the number, the better.



We also like to recommend practicing on a piece of chicken or steak, so you get a feel for the power levels.

We do not sell smoke evacuators for this. The amount of smoke is very little.

If you have any more questions, please reach out to your local EICKEMEYER representative.