

The Best Intensive Care for your patients

VETARIO T40M & T50M



Give your patients the best possible chance of recovery with Vetario Intensive Care Units. Trusted by veterinarians, breeders, and wildlife rescuers, our enclosures are designed to provide a **safe, warm, and hygienic environment** for critically ill animals and birds.

Designed for optimal patient care, the T40M and T50M units offer comprehensive support with integrated temperature control, humidity management, drug nebulisation, and oxygen therapy.

They are robust, reliable, easy to use and can be quickly set up for emergency cases.

KEY BENEFITS FOR VETERINARIANS AND THEIR PATIENTS



ENHANCED RECOVERY WITH OXYGEN THERAPY

These units feature an oxygen connector for safe operation at elevated oxygen levels, which can **significantly accelerate patient recovery during critical care.**



ENHANCED PATIENT MONITORING

Tough, clear door with a secure latch ensures patient safety while allowing easy observation. Switchable light to **minimise patient disturbance during monitoring.**



MINIMISING PATIENT RECOVERY TIME

Using a Vetario ICU enables more efficient patient warming, therefore, **speeding up recovery – saving your busy surgery time and money!**



SAFETY AND PEACE OF MIND

Advanced control system: includes visual and audible alerts for internal and ambient temperature deviations, **ensuring patient safety.**



COMPREHENSIVE ENVIRONMENTAL CONTROL

– digital display showing real-time temperature and humidity levels for **enhanced comfort and recovery.**



EFFORTLESS DRUG ADMINISTRATION


Deliver medication conveniently through the supplied medicine reservoir, compatible with nebuliser pumps available from your Vetario distributor, streamlining treatment protocols.




HYGIENIC & EASY TO CLEAN:

smooth interior with curved edges and an option to disassemble, **making deep cleaning easy and efficient.**

OXYGEN SAFE DESIGN

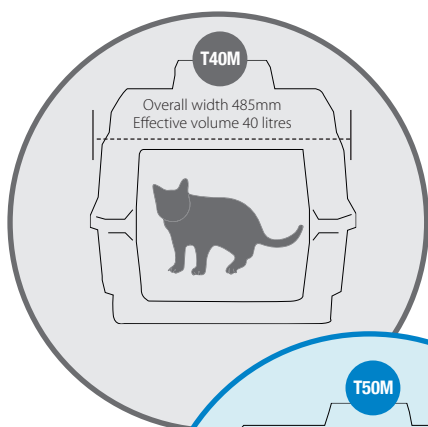

The T40M is suitable for rodents, kittens and small puppies, etc. The larger T50M can accommodate a small dog or queen and her litter, for example.

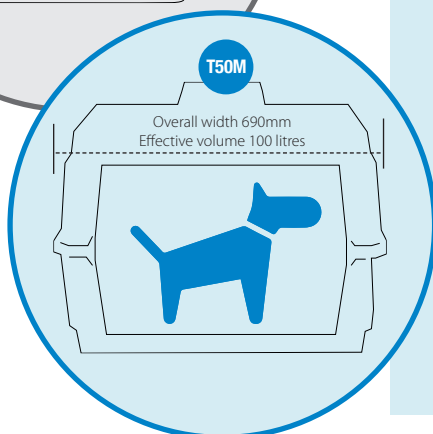



'Vetario units quickly became essential across all areas of our practice, meeting the needs of dogs, cats, and exotic animals alike. They are especially valuable in recovery areas, where maintaining oxygen levels and regulating body temperature are critical. In our intensive care room, Vetario units provide constant monitoring and support for our most critical patients, including those recovering from anaemia, trauma, or respiratory distress' **Centre Hospitalier Vétérinaire, Saint-Martin**




Why choose Vetario?

-  Accurate and adjustable temperature control
-  Oxygen compatibility
-  Smooth, easy-to-clean interior
-  No vulnerable edges for animals or birds to chew
-  Tough, clear hinged door with secure catch
-  Positive pressure with gentle fan ventilation
-  Twin stage filter, for removal of harmful bacteria and fungi
-  Water reservoir for increasing humidity with external top up
-  Optional EX Humidity Pump integration enables fully automated humidity control
-  Digital temperature control
-  Nebuliser ampule fitting
-  Switchable light
-  Low weight makes them easy to carry
-  Cabinet protected by Biomaster Product Protection

Oxygen Safety in Veterinary ICUs

Elevated oxygen levels are crucial for animal recovery but increase fire risk when combined with electronics, fans, or heaters. These units eliminate all potential ignition sources in oxygen-enriched environments, even in single-fault conditions such as leaks or component failures. The T40M and T50M models are rigorously tested by the British Standards Institute (BSI) and meet stringent safety standards, classified as RX products under USA FDA guidelines. This ensures veterinarians can safely utilize oxygen therapy without compromising patient or clinic safety.

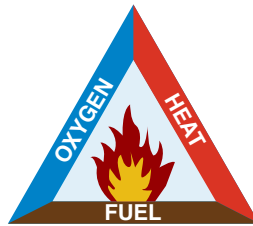

SPECIFICATIONS	T40M	T50M
Oxygen compatible	•	•
Digital temperature control	•	•
Display choice of °C or °F	•	•
Digital Humidity display	•	•
Water pan for humidity	•	•
EX Humidity pump option	•	•
Temperature alarm	•	•
Room temperature alarm	•	•
Secondary thermostat	•	•
LED light	•	•
Nebuliser fitting	•	•
Antimicrobial plastic	•	•
Overall height (mm)	470	550
Overall width (mm)	485	690
Overall depth (mm)	385	490
Floor area (mm)	400 x 300	600 x 400
Effective volume (litres)	40	100
Power consumption (Watts) Max	150	150
Typical	85	100
Power supply	115v	115v
Weight (kg)	6.7	8.7

All filters are washable and replacement filter/evaporating block sets are available from Brinsea or your local stockist.

In order that we can continue our policy of innovation and improvement we reserve the right to alter specifications without notice. Vetario intensive care units are subject to EU Design Registration No. 001779729-0005.

'Vetario incubators are exceptional for patient recovery, halving recovery time by providing a safe and controlled environment. This not only saves time but also increases efficiency, allowing more patients to recover safely and cost-effectively.'

Sonya Miles, BVSc CertAVP(ZM), CertAqV, MRCVS. RCVS Recognised Advanced Practitioner in Zoological Medicine, WAVMA Certified Aquatic Species Veterinarian



'I would not be without the Vetario units and they have certainly helped in saving many lives. The practicality makes life much easier and eliminates the concerns I would have had previously regarding thermoregulation.'

Jessie Blyth RVN, APVN (Wildlife)



BRINSEA PRODUCTS LTD
32-33 Buckingham Road,
Weston Industrial Estate,
Weston-super-Mare,
BS24 9BG
Tel: +44 (0) 1934 417523

GUARANTEE: A comprehensive **FREE 3 year extended guarantee** is available on all Vetario branded intensive care units when you register online within 30 days. Should a fault arise during the guarantee period the unit may be repaired or replaced at Brinsea's discretion.
www.vetario.co.uk or www.vetario.com (USA & Canada)

