




EICKEMEYER® NarkoVet
Anaesthetic Unit

veterinary technology for life
Eickemeyer®

More flexibility during anaesthesia

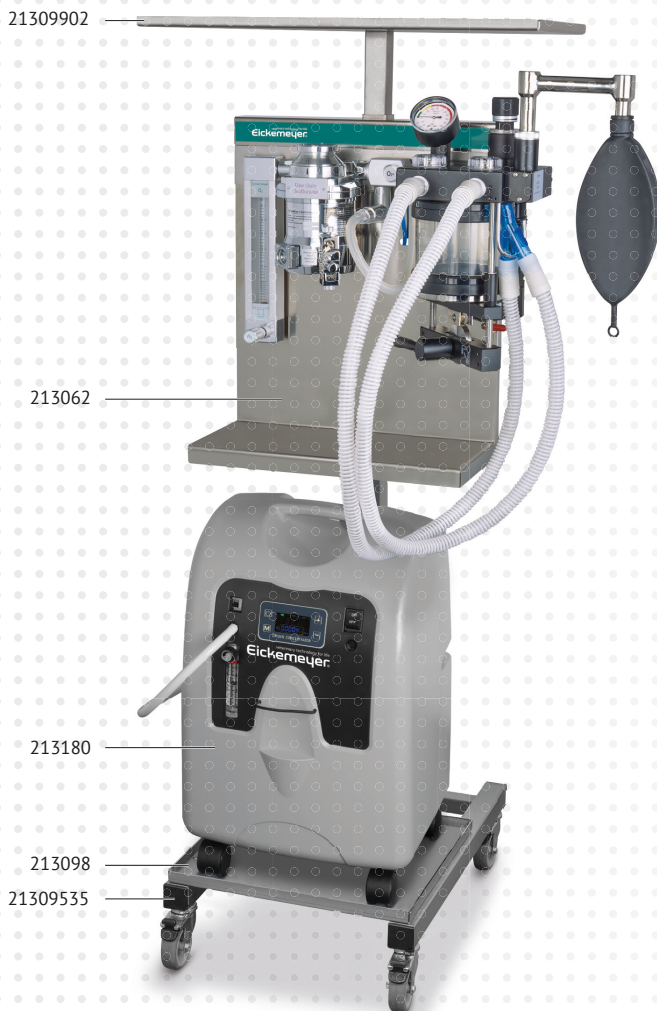

EICKEMEYER® NarkoVet Anaesthetic Unit

The EICKEMEYER® NarkoVet Anaesthetic Unit is now even more flexible: Can be set up as a desk, wall or stand model.

Table Mounted Model

- ▶ Flow Meter (0.1 – 5.0 l/min)
- ▶ System pressure gauge
- ▶ Infinitely adjustable pressure relief valve
- ▶ Oxygen bypass
- ▶ TEC 3 Isoflourane Vaporiser
- ▶ Soda Lime canister and accessories
- ▶ Dimensions (in cm): L 32 x W 25 x H 51
- ▶ Suitable for small animals up to 80 kg
- ▶ OxyVet IV Oxygen Concentrator (213180) sold separately

213062

Wall Bracket for EICKEMEYER® Anaesthetic Unit

- ▶ Enables the EICKEMEYER® Anaesthetic Unit to be mounted to the wall

21309403

Anaesthetic Unit Trolley

21309535

Anaesthetic Unit Oxygen Concentrator Shelf

213098

Anaesthetic Unit Top Shelf for monitoring devices

- ▶ Stainless steel
- ▶ Dimensions: 60 x 30 cm

21309902

EICKEMEYER®

3 Windmill Business Village
Brooklands Close, Sunbury-on-Thames
Surrey, TW16 7DY, United Kingdom

T 020 8891 2007
info@eickemeyer.co.uk

TELEPHONE 020 8891 2007

www.eickemeyer.co.uk



EICKEMEYER® NarkoVet Anaesthetic Unit

veterinary technology for life

Eickemeyer®



213180

NarkoVet Accessories

OxyVet IV Oxygen concentrator

Oxygen concentrator for use in oxygen therapy and anesthesia. An oxygen concentrator draws in room air, separates the oxygen from other gases in the environment and delivers the oxygen at high concentrations.

- ▶ Saving costs and effort: No oxygen bottle required, lower power consumption
- ▶ Simple operation thanks to the clear arrangement of the setting functions and monitoring information
- ▶ High patient safety through visual and acoustic alarms
- ▶ Noise level: ≤ 55 dbA
- ▶ Flow rate continuously adjustable from 0.5 – 5 l/min
- ▶ Output of up to 93 % oxygen (± 3 %)
- ▶ Incl. Bubble humidifier and O₂ hose
- ▶ Dimensions (in cm): L 34 x W 37 x H 67
- ▶ Weight: 15.3 kg
- ▶ 1 year warranty

213180

Isoflurane Safety Key Filler

213246

EickSorber Anaesthetic Gas Absorber

- ▶ Includes 2 m exhaust hose

213488

Soda Lime Anaesthetic Absorber

- ▶ Low dust particles
- ▶ Colour indicator: white to violet
- ▶ 5 Litre

215206



213488



215206

EICKEMEYER®

3 Windmill Business Village
Brooklands Close, Sunbury-on-Thames
Surrey, TW16 7DY, United Kingdom

T 020 8891 2007
F 020 8891 2686
info@eickemeyer.co.uk

TELEPHONE 020 8891 2007

www.eickemeyer.co.uk